


LABORATORY EQUIPMENT

Product catalogue

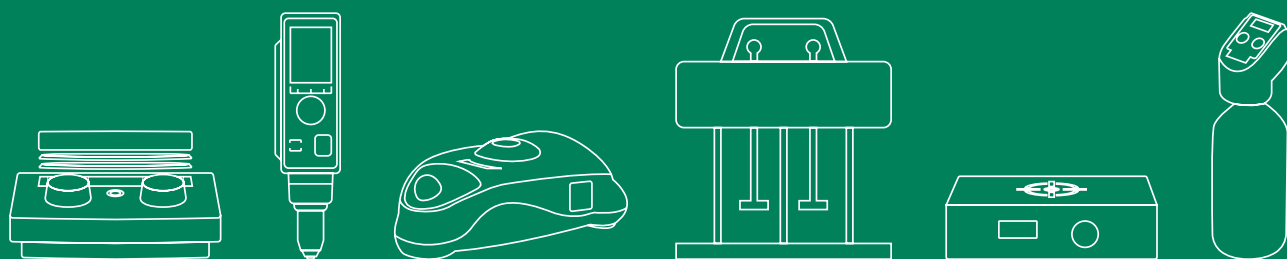

VELP Scientifica

Established in 1983, VELP Scientifica is today one of the world's leading manufacturers of analytical instruments and laboratory equipment that has made an impact on the world-wide market with Italian products renowned for innovation, design and premium connectivity.

We owe our success to a rich patrimony of ideas, significant investments in R&D and technological development aiming at maximizing productivity as well as reducing time-to-market.

VELP is a pioneer of smart lab technology. Our analytical instruments and laboratory equipment are accelerating scientific work every day giving scientists the possibility to monitor and control VELP technology from a virtual dashboard accessible from any device, everywhere and anytime!

VELP is today a global company with an expanding local focus. Our internationalization process has led us to open local subsidiaries that provide service to our customers in different parts of the world:

- **VELP Scientific, Inc.** in Long Island (USA) offering sales and technical support for the United States, Canada, and Mexico.
- **VELP China Co. LTD** in Shanghai guarantees a local presence in China through commercial and technical support.

VELP has always generated value for its customers through close cooperation with its qualified distribution network which today counts more than 300 distributors in 130 countries.

Product Overview

VELP Ermes

- SMART LAB CLOUD PLATFORM 7
- VELP ERMES SUBSCRIPTIONS 7

Magnetic Stirrers

- MSL 50 Digital 8
- MSL 25 Digital 8
- MSL 8 Digital 9
- MSL 8 9
- MST Digital 10
- MST 10
- MICROSTIRRER 11
- ESP 11
- MULTISTIRRER 15 Digital 12
- MULTISTIRRER 6 Digital 12
- MULTISTIRRER 15 13
- MULTISTIRRER 6 13
- AMI 14
- AMI 4 14

Heating Magnetic Stirrers

- AREX-6 Connect PRO 16
- AREX-6 Digital PRO 17
- AREX-6 18
- ARE-6 19
- AREX Digital PRO 20
- AREX 21
- ARE 22
- AM4 Digital PRO 23
- AM4 24
- AREC Connect 27
- AREC.X 28
- AREC 29
- HSC 29
- MULTI-HS 15 Digital 31
- MULTI-HS 6 Digital 32
- VTF 33

Heating Plates

| | |
|-------|----|
| ▪ RC | 34 |
| ▪ RC2 | 34 |
| ▪ REC | 35 |

Overhead Stirrers

| | |
|----------------------|----|
| ▪ OHS 200 Advance | 37 |
| ▪ OHS 100 Advance | 37 |
| ▪ OHS 60 Advance | 38 |
| ▪ OHS 200 Digital | 38 |
| ▪ OHS 100 Digital | 39 |
| ▪ OHS 60 Digital | 39 |
| ▪ OHS 40 Digital | 40 |
| ▪ OHS 20 Digital | 40 |
| ▪ DLH | 41 |
| ▪ DLS | 41 |
| ▪ PW | 42 |
| ▪ LH | 42 |
| ▪ LS | 43 |
| ▪ ES | 43 |
| ▪ CONTROLLER Advance | 46 |
| ▪ CONTROLLER Digital | 46 |

Vortex Mixers

| | |
|-------------|----|
| ▪ Wizard IR | 47 |
| ▪ Classic | 47 |
| ▪ TX4 | 48 |
| ▪ ZX4 | 48 |
| ▪ ZX3 | 49 |
| ▪ RX3 | 49 |
| ▪ MULTI-TX5 | 51 |

Homogenizer

| | |
|-------|----|
| ▪ OV5 | 52 |
|-------|----|

COD Thermoreactors

| | |
|----------|----|
| ▪ ECO 16 | 54 |
| ▪ ECO 6 | 54 |
| ▪ ECO 25 | 55 |
| ▪ ECO 8 | 55 |

BOD and Respirometers

| | |
|--|----|
| ▪ RESPIROMETRIC Sensor System 6 BOD | 57 |
| ▪ RESPIROMETRIC Sensor System 6 Maxi BMP | 57 |
| ▪ RESPIROMETRIC Sensor System 6 for Soil Analysis | 58 |
| ▪ RESPIROMETRIC Sensor System 6 for Plastic Biodegradability | 58 |
| ▪ BOD Sensor | 59 |
| ▪ BOD Sensor System 6 | 59 |
| ▪ BOD Sensor System 10 | 60 |

Cooled Incubators

| | |
|---------------------|----|
| ▪ FOC 200IL Connect | 61 |
| ▪ FOC 200I Connect | 61 |
| ▪ FOC 200E Connect | 62 |
| ▪ FOC 120I Connect | 62 |
| ▪ FOC 120E Connect | 63 |
| ▪ FTC 120 | 63 |

Flocculators

| | |
|---------|----|
| ▪ JLT 6 | 65 |
| ▪ JLT 4 | 65 |
| ▪ FC6S | 66 |
| ▪ FC4S | 66 |
| ▪ FP4 | 67 |


Other Solutions

| | |
|---------------------------------------|----|
| ▪ ROTAX 6.8 Overhead Mixer | 68 |
| ▪ OCB Open Circulating Baths | 68 |
| ▪ Radiation Detector | 69 |
| ▪ SP 311 Peristaltic Pumps | 69 |
| ▪ TB1 Turbidimeter | 70 |
| ▪ JPV Recirculating Water Vacuum Pump | 70 |

VELP Ermes

■ SMART LAB CLOUD PLATFORM

VELP Ermes creates a connected ecosystem of devices, people and data that transmit information between each other cutting down distances and expanding your scientific potential.

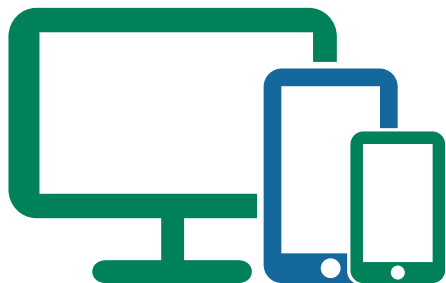

Monitor and Control your instruments

Monitor and manage multiple instruments 24/7 from your internet devices. Real-time visualization of your analysis and of the instrument working conditions. Immediate notifications and alerts.



Remote Service Support

thanks to the direct contact with VELP specialists and official partners. With VELP Ermes we can guide you to get the maximum from your instrument.



Safely Manage your Data

Collect and store your data with the maximum level of encryption and cyber-security standards. Generate and share reports of your analysis in multiple formats.



Your Instrument Always Updated

Effortlessly connect your instrument via secure Wi-Fi or cable and update the software immediately and free of charge.


■ VELP ERMES SUBSCRIPTIONS

Configure your account and start enjoying all the services and capabilities of VELP Ermes and your instrument for free for 3 months. Decide later whether to extend your subscription with an Ermes Card for 1 and 3 years. Each Ermes Card allows you to connect up to 10 instruments in your company.




3
MONTHS
FREE TRIAL

FREE TRIAL
3 MONTHS PER INSTRUMENT

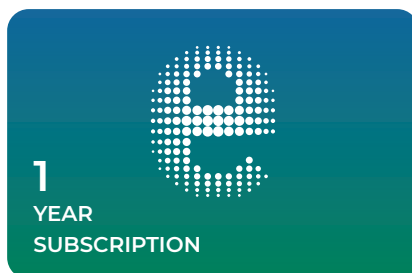

1
YEAR
SUBSCRIPTION

1 YEAR / 10 INSTRUMENTS
Code E00010012

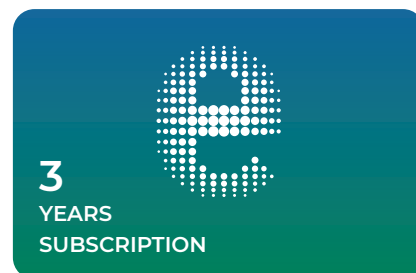

3
YEARS
SUBSCRIPTION

3 YEARS / 10 INSTRUMENTS
Code E00010036

Magnetic Stirrers

■ MSL 50 Digital

High volume magnetic stirrer with digital display to accurately set the working conditions. Programmable operations with the integrated timer and auto-reverse. Main applications of MSL 50 Digital are high volumes preparations, cell cultures and tissue cultures, buffer solution preparations.

PRECISE AND PROGRAMMABLE

SpeedServo™ technology maintains constant speed. Unattended stirring from 1 minute to 99 h 59 minutes.

ROBUST AND VERSATILE

The brushless motor provides continuous operation for high volumes at any speed. Robust stainless steel top for easy visualization of color changes and particles.



Volume (H₂O)
50 L



Stirring speed
From 20 to 1000 rpm



Timer &
Auto-reverse



Safety Lock



INSTRUMENT - CODE

MSL 50 Digital

100÷240 V / 50-60 Hz

F203A0512



■ MSL 25 Digital

High volume magnetic stirrer with digital display to accurately set the working conditions. Programmable operations with the integrated timer and auto-reverse. Main applications of MSL 25 Digital are high volumes preparations, cell cultures and tissue cultures, buffer solution preparations.

PRECISE AND PROGRAMMABLE

SpeedServo™ technology maintains constant speed. Unattended stirring from 1 minute to 99 h 59 minutes.

ROBUST AND VERSATILE

The brushless motor provides continuous operation for high volumes at any speed. Robust stainless steel top for easy visualization of color changes and particles.



Volume (H₂O)
25 L



Stirring speed
From 40 to 1500 rpm



Timer &
Auto-reverse



Safety Lock



INSTRUMENT - CODE

MSL 25 Digital

100÷240 V / 50-60 Hz

F203A0511



■ MSL 8 Digital

High volume magnetic stirrer with digital display to accurately set the working conditions. Programmable operations with the integrated timer and auto-reverse. Main applications of MSL 8 Digital are high volumes preparations, cell cultures and tissue cultures, buffer solution preparations.

PRECISE AND PROGRAMMABLE

SpeedServo™ technology maintains constant speed. Unattended stirring from 1 minute to 99 h 59 minutes.

ROBUST AND VERSATILE

The brushless motor provides continuous operation for high volumes at any speed. Robust stainless steel top for easy visualization of color changes and particles.



Volume (H₂O)
8 L



Stirring speed
From 40 to 1500 rpm



Timer &
Auto-reverse



Safety Lock



INSTRUMENT - CODE

MSL 8 Digital

100÷240 V / 50-60 Hz

F203A0510



■ MSL 8

Magnetic stirrer guaranteeing continuous operations at any speed ensuring safety and resistance to chemical agents. MSL 8 has a gentle start in order to protect the sample. The control panel is protected from liquid spills.

PRECISE

SpeedServo™ technology maintains constant speed.

ROBUST AND VERSATILE

The brushless motor provides continuous operation for high volumes at any speed. Robust stainless steel top for easy visualization of color changes and particles.



Volume (H₂O)
8 L



Stirring speed
From 40 to 1500 rpm



INSTRUMENT - CODE

MSL 8

100÷240 V / 50-60 Hz

F203A0513



MST Digital

Compact magnetic stirrer with digital display with accurate speed setting, timer and auto-reverse. The white surface is designed to facilitate the visualization of particles and color change. MST Digital allows continuous operation and remains cold even after days of continuous use. Main applications in microbiology, biochemistry and chemistry.

PRECISE AND PROGRAMMABLE

SpeedServo™ technology maintains constant speed. Unattended stirring from 1 to 900 minutes.

POWERFUL

Equipped with a powerful brushless motor.



Volume (H₂O)
5 L



Stirring speed
Up to 1500 rpm



Timer &
Auto-reverse



Safety Lock



INSTRUMENT - CODE

MST Digital

100÷240 V / 50-60 Hz

F203A0450



MST

Powerful and compact magnetic stirrer with a low profile. The white surface is designed to facilitate the visualization of particles and color change. MST Digital allows continuous operation and remains cold even after days of continuous use. Main applications in microbiology, biochemistry and chemistry.

RESISTANT

The technopolymer structure ensures maximum resistance to chemicals.



Volume (H₂O)
5 L



Stirring speed
Up to 1100 rpm



INSTRUMENT - CODE

MST

100÷240 V / 50-60 Hz

F203A0440



■ MICROSTIRRER

Compact and light weight magnetic stirrer with metal structure for visualization of particles and color changes. Specially designed for microtitrations.

RELIABLE

It remains cold, even after days of continuous use and guarantees excellent speed control, even at low speed.



Volume (H₂O)
5 L



Stirring speed
Up to 1100 rpm



INSTRUMENT - CODE

MICROSTIRRER

100÷240 V / 50-60 Hz

F203A0161



■ ESP

Flat and light weight magnetic stirrer with no moving mechanical components: it is therefore maintenance-free. The clear surface is ideal for visualization of particles and color changes.

RESISTANT

The control panel is protected from liquid spills with run-off groove and the technopolymer casing resists chemicals.



Volume (H₂O)
5 L



Stirring speed
Up to 1100 rpm



INSTRUMENT - CODE

ESP

100÷240 V / 50-60 Hz

F206A0179



MULTISTIRRER 15 Digital

Multi position magnetic stirrer with stainless steel top for the simultaneous stirring of 15 positions. A bright digital display to accurately set the stirring speed, timer, and auto-reverse. Ideal for chemistry, microbiology, biochemistry applications.

PRECISE AND PROGRAMMABLE

SpeedServo™ technology maintains constant speed. Unattended stirring from 1 to 900 minutes.

RESISTANT

Stainless steel top plate. The control panel is shielded from liquid spills.



Volume (H₂O)
250 ml (x15)



Stirring speed
From 80 to 1500 rpm



Timer &
Auto-reverse



Safety Lock



INSTRUMENT - CODE

MULTISTIRRER 15 Digital

100÷240 V / 50-60 Hz

F203A0180



MULTISTIRRER 6 Digital

Multi position magnetic stirrer with stainless steel top for the simultaneous stirring of 6 positions. A bright digital display to accurately set the stirring speed, timer, and auto-reverse. Ideal for chemistry, microbiology, biochemistry applications.

PRECISE AND PROGRAMMABLE

SpeedServo™ technology maintains constant speed. Unattended stirring from 1 to 900 minutes.

RESISTANT

Stainless steel top plate. The control panel is shielded from liquid spills.



Volume (H₂O)
400 ml (x6)



Stirring speed
From 80 to 1500 rpm



Timer &
Auto-reverse



Safety Lock



INSTRUMENT - CODE

MULTISTIRRER 6 Digital

100÷240 V / 50-60 Hz

F203A0179



MULTISTIRRER 15

Multi position magnetic stirrer with stainless steel top for the simultaneous stirring of 15 positions. Pale surface for visualization of particles and color changes. Remains cold, even after days of continuous use. Ideal for chemistry, microbiology, biochemistry applications.

PRECISE

SpeedServo™ technology maintains constant speed.

RESISTANT

Stainless steel top plate. The control panel is shielded from liquid spills.



Volume (H₂O)
250 ml (x15)



Stirring speed
From 80 to 1500 rpm



INSTRUMENT - CODE

MULTISTIRRER 15

100÷240 V / 50-60 Hz

F203A0178



MULTISTIRRER 6

Multi position magnetic stirrer with stainless steel top for the simultaneous stirring of 6 positions. Pale surface for visualization of particles and color changes. Remains cold, even after days of continuous use. Ideal for chemistry, microbiology, biochemistry applications.

PRECISE

SpeedServo™ technology maintains constant speed.

RESISTANT

Stainless steel top plate. The control panel is shielded from liquid spills.



Volume (H₂O)
400 ml (x6)



Stirring speed
From 80 to 1500 rpm



INSTRUMENT - CODE

MULTISTIRRER 6

100÷240 V / 50-60 Hz

F203A0177



■ AMI

Illuminated magnetic stirrer particularly useful for titrations. The illuminated surface facilitates the visualization of particles and color changes. It guarantees excellent speed control (electronic), even at low speed.

RESISTANT

Epoxy-coated tough steel casing and stainless steel cover. The control panel is shielded from liquid spills

RELIABLE

It remains cold, even after days of continuous use and is equipped with a pale surface for visualization of particles and color changes.



Volume (H₂O)
5 L



Stirring speed
Up to 1100 rpm



LED



INSTRUMENT - CODE

AMI

100÷240 V / 50-60 Hz

F204A0167



■ AMI 4

Illuminated multi-position magnetic stirrer particularly useful for titrations. The illuminated surface facilitates the visualization of particles and color changes. It guarantees excellent speed control (electronic), even at low speed.

RESISTANT

Epoxy-coated tough steel casing and stainless steel cover. The control panel is shielded from liquid spills

RELIABLE

It remains cold, even after days of continuous use and is equipped with a pale surface for visualization of particles and color changes.



Volume (H₂O)
5 L (x4)



Stirring speed
Up to 1100 rpm



LED



INSTRUMENT - CODE

AMI 4

100÷240 V / 50-60 Hz

F204A0168



ACCESSORIES - MAGNETIC STIRRERS

| | MST | MST Digital | MICRO-STIRRER | ESP | MSL 8 | MSL 8 Digital | MSL 25 Digital | MSL 50 Digital | AMI | AMI 4 | MULTI-STIRRER | MULTI-STIRRER Digital |
|--|-----|-------------|---------------|-----|-------|---------------|----------------|----------------|-----|-------|---------------|-----------------------|
| A00001057 Magnetic stir bar, Ø 6x20 mm | • | • | • | | | | | | • | • | • | • |
| A00001056 Magnetic stir bar, Ø 6x35 mm | • | • | • | | • | • | | | • | • | • | • |
| A00001061 Magnetic stir bar, Ø 9.5x60 mm | | | | | | | | • | | | | |
| A00000356 Magnetic stir bar, Ø 8x40 mm | | | | • | • | • | | | | | | |
| A00000387 Magnetic stir bar, Ø 8x50 mm | | | | | | | • | | | | | |
| A00001055 Thermostatic bath for samples, 280x180x410 mm | | | | | | | | | | | • | • |
| 10003083 Plug USA for power supply | • | • | • | • | • | • | • | • | • | • | • | • |
| 10003084 Plug UK for power supply | • | • | • | • | • | • | • | • | • | • | • | • |
| 10003085 Plug AU for power supply | • | • | • | • | • | • | • | • | • | • | • | • |



Heating Magnetic Stirrers

■ AREX-6 Connect PRO

Digital hot plate stirrer that connects effortlessly to the VELP Ermes cloud platform via secure Wi-Fi allowing monitoring of and control of the reactions at any time. The white LED display and a set of icons provide precise settings and clear readout of the working conditions. The fast heating, 135 mm diameter CerAlTop™ plate protects from chemicals, scratches and enhances the lifetime of the unit.

SAFETY AND DURABILITY

Two independent safety circuits protect the instrument in case of over-temperature, switching-off the heating. IP42 certified.

STATE-OF-THE-ART THERMOREGULATION

Extremely stable and fast thermoregulation with minimal overshoot thanks to the PID thermoregulation software.

SUPERIOR PERFORMANCES

Brushless motor and Alnico Magnet provide state-of-the-art stirring performance.



ermes enabled



Volume (H₂O)
20 L



Stirring speed
From 30 to 1700 rpm



Temperature range
From room temperature
to 370 °C



Timer &
Auto-reverse




Safety Lock



External connection
Pt 100 Probe

CONFIGURATIONS - AREX-6 Connect PRO

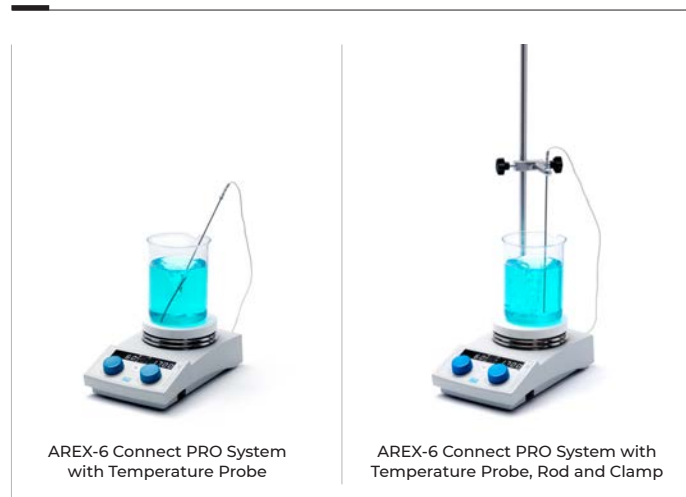

AREX-6 Connect PRO System
with Temperature Probe

AREX-6 Connect PRO System with
Temperature Probe, Rod and Clamp

INSTRUMENT - CODE

| | | |
|--|------------------|------------|
| AREX-6 Connect PRO with Temperature Probe | 230 V / 50-60 Hz | SA20500465 |
| AREX-6 Connect PRO with Temperature Probe | 115 V / 60 Hz | SA20510465 |
| AREX-6 Connect PRO with Temperature Probe, Rod and Clamp | 230 V / 50-60 Hz | SB20500465 |
| AREX-6 Connect PRO with Temperature Probe, Rod and Clamp | 115 V / 60 Hz | SB20510465 |



■ AREX-6 Digital PRO

Digital hot plate stirrer with white LED display and a set of icons providing precise settings and clear readout of the working conditions. The fast heating, 135 mm diameter CerAlTop™ plate protects from chemicals, scratches and enhances the lifetime of the unit. Enhanced kinetics and programmable operations thanks to integrated timer and auto-reverse.

SAFETY AND DURABILITY

Two independent safety circuits protect the instrument in case of over-temperature, switching-off the heating. IP42 certified.

STATE-OF-THE-ART THERMOREGULATION

Extremely stable and fast thermoregulation with minimal overshoot thanks to the PID thermoregulation software.

SUPERIOR PERFORMANCES

Brushless motor and Alnico Magnet provide state-of-the-art stirring performance.



Volume (H₂O)
20 L



Stirring speed
From 30 to 1700 rpm



Temperature range
From room temperature
to 370 °C



Timer &
Auto-reverse



Safety Lock



External connection
Pt 100 Probe



External connection
VTF

CONFIGURATIONS - AREX-6 Digital PRO



AREX-6 Digital PRO System
with Temperature Probe

AREX-6 Digital PRO System with
Temperature Probe, Rod and clamp

AREX-6 Digital PRO System with
VTF Thermoregulator

INSTRUMENT - CODE

| | | |
|--|------------------|------------|
| AREX-6 Digital PRO with Temperature Probe | 230 V / 50-60 Hz | SA20500460 |
| AREX-6 Digital PRO with Temperature Probe | 115 V / 60 Hz | SA20510460 |
| AREX-6 Digital PRO with Temperature Probe, Rod and Clamp | 230 V / 50-60 Hz | SD20500460 |
| AREX-6 Digital PRO with Temperature Probe, Rod and Clamp | 115 V / 60 Hz | SD20510460 |
| AREX-6 Digital PRO with VTF | 230 V / 50-60 Hz | SB20500460 |
| AREX-6 Digital PRO with VTF | 115 V / 60 Hz | SB20510460 |



■ AREX-6

Safe and powerful hot plate stirrer with 135 mm diameter CerAlTop™ plate that protects from chemicals, scratches and enhances the lifetime of the unit. Combinable with VTF for accurate thermoregulation of the sample ± 0.5 °C.

SAFETY AND DURABILITY

The illuminated heating switch on/off prevents unintentional heat up and the illuminated heat icon informs the operator when the top plate temperature is above 50°C.

SUPERIOR PERFORMANCES

Brushless motor and Alnico Magnet provide state-of-the-art stirring performance.



Volume (H₂O)
20 L



Stirring speed
Up to 1100 rpm



Temperature range
From room temperature
to 370 °C



External connection
VTF



CONFIGURATIONS - AREX-6



AREX-6

AREX-6 System with
VTF Thermoregulator

INSTRUMENT - CODE

| | | |
|-----------------|------------------|------------|
| AREX-6 | 230 V / 50-60 Hz | F20500463 |
| AREX-6 | 115 V / 60 Hz | F20510463 |
| AREX-6 with VTF | 230 V / 50-60 Hz | SA20500463 |
| AREX-6 with VTF | 115 V / 60 Hz | SA20510463 |



■ ARE-6

Safe and powerful hot plate stirrer with 135 mm diameter aluminium plate. Precise analog knobs to control speed and plate temperature with the maximum ease-of-use.

SAFETY AND DURABILITY

The illuminated heating switch on/off prevents unintentional heat up and the illuminated heat icon informs the operator when the top plate temperature is above 50 °C.

SUPERIOR PERFORMANCES

Brushless motor and Alnico Magnet provide state-of-the-art stirring performance.



Volume (H₂O)
20 L



Stirring speed
From 30 to 1700 rpm



Temperature range
From room temperature
to 370 °C



INSTRUMENT - CODE

| | | |
|-------|------------------|-----------|
| ARE-6 | 230 V / 50-60 Hz | F20500462 |
| ARE-6 | 115 V / 60 Hz | F20510462 |

ACCESSORIES - ARE/AREX-6 Series



| | AREX-6 CONNECT PRO | AREX-6 Digital PRO | AREX-6 | ARE-6 |
|---|--------------------|--------------------|--------|-------|
| A00001062 Magnetic stir bar, Ø 3x6 mm | • | • | • | • |
| A00001063 Magnetic stir bar, Ø 4,5x12 mm | • | • | • | • |
| A00001057 Magnetic stir bar, Ø 6x20 mm | • | • | • | • |
| A00001056 Magnetic stir bar, Ø 6x35 mm | • | • | • | • |
| A00000356 Magnetic stir bar, Ø 8x40 mm | • | • | • | • |
| A00001061 Magnetic stir bar, Ø 10x60 mm | • | • | • | • |
| A00000336 Magnetic cross shape stir bar, Ø 10x5 mm | • | • | • | • |
| A00000352 Magnetic cross shape stir bar, Ø 20x8 mm | • | • | • | • |
| A00000354 Magnetic disc stir bar, Ø 10x6 mm | • | • | • | • |
| A00000355 Magnetic disc stir bar, Ø 20x10 mm | • | • | • | • |
| A00000357 Magnetic octagonal stir bar, Ø 10x13 mm | • | • | • | • |
| A00000268 Temperature Probe AISI 316 Ti Ø 3 mm | • | • | • | • |
| A00000349 Temperature probe glass-coated Ø 5 mm | • | • | • | • |
| F208B0063 VTF Digital Thermoregulator | | • | • | |
| A00000382 Extension for support rod | • | • | • | • |
| A00000280 Clamp for the probe | • | • | • | • |
| A00000351 Handle for AluBlock removal | • | • | • | • |
| A00001069 Support rod | • | • | • | • |
| E00010012 VELP Ermes 1 years Connection | • | | | |
| E00010036 VELP Ermes 3 years Connection | • | | | |

■ AREX Digital PRO

Digital hot plate stirrer ensuring precise settings and clear readout of the working conditions. The fast heating, 135 mm diameter CerAlTop™ plate protects from chemicals, scratches and enhances the lifetime of the unit.

SAFETY AND DURABILITY

The control panel is protected against possible damage caused by liquid spills thanks to a dedicated run-off groove. IP 42 certified.

PRECISE THERMOREGULATION

Excellent thermoregulation performance with Pt100 Probe and VTF.



Volume (H₂O)
20 L



Stirring speed
Up to 1500 rpm



Temperature range
From room temperature
to 370 °C



Timer



External connection
Pt 100 Probe



External connection
VTF

CONFIGURATIONS - AREX Digital PRO



AREX Digital PRO System with
Temperature Probe

AREX Digital PRO System with
Temperature Probe, Rod and clamp

AREX Digital PRO System with
VTF Thermoregulator

INSTRUMENT - CODE

| | | |
|--|------------------|------------|
| AREX Digital PRO with Temperature Probe | 230 V / 50-60 Hz | SA20500410 |
| AREX Digital PRO with Temperature Probe | 115 V / 60 Hz | SA20510410 |
| AREX Digital PRO with Temperature Probe, Rod and Clamp | 230 V / 50-60 Hz | SC20500410 |
| AREX Digital PRO with Temperature Probe, Rod and Clamp | 115 V / 60 Hz | SC20510410 |
| AREX Digital PRO with VTF | 230 V / 50-60 Hz | SB20500410 |
| AREX Digital PRO with VTF | 115 V / 60 Hz | SB20510410 |



■ AREX

Reliable hot plate stirrer with 135 mm diameter CerAlTop™ plate that protects from chemicals, scratches and enhances the lifetime of the unit. Combinable with VTF for accurate thermoregulation of the sample ± 0.5 °C.

SAFETY AND DURABILITY

The control panel is protected against possible damage caused by liquid spills thanks to a dedicated run-off groove. IP 42 certified.

PRECISE THERMOREGULATION

Excellent thermoregulation performance with VTF.



Volume (H₂O)
20 L



Stirring speed
Up to 1500 rpm



Temperature range
From room temperature
to 370 °C



External connection
VTF



CONFIGURATIONS - AREX



AREX

AREX System with
VTF Thermoregulator

INSTRUMENT - CODE

| | | |
|---------------|------------------|------------|
| AREX | 230 V / 50-60 Hz | F20500413 |
| AREX | 115 V / 60 Hz | F20510413 |
| AREX with VTF | 230 V / 50-60 Hz | SA20500413 |
| AREX with VTF | 115 V / 60 Hz | SA20510413 |



■ ARE

Hot plate stirrer with SpeedServo™ that ensures constant speed even when the viscosity changes. ARE guarantees an excellent speed control and it is widely used in research and development, industrial and university laboratories worldwide.

DURABLE

The ARE is IP42 certified and thanks to the resistant die-cast housing and the run-off groove all the internal parts of the instrument are protected from potential damages.

RELIABLE

Aluminum alloy top plate (135 mm diameter) ensures a great temperature homogeneity and optimum heat transfer across the entire surface.



Volume (H₂O)
15 L



Stirring speed
Up to 1500 rpm



Temperature range
From room temperature
to 370 °C

INSTRUMENT - CODE

| | | |
|-----|------------------|-----------|
| ARE | 230 V / 50-60 Hz | F20500162 |
| ARE | 115 V / 60 Hz | F20510162 |



ACCESSORIES - ARE-AREX Series

| | AREX Digital PRO | AREX | ARE |
|---|------------------|------|-----|
| A00001062 Magnetic stir bar, Ø 3x6 mm | • | • | • |
| A00001063 Magnetic stir bar, Ø 4,5x12 mm | • | • | • |
| A00001057 Magnetic stir bar, Ø 6x20 mm | • | • | • |
| A00001056 Magnetic stir bar, Ø 6x35 mm | • | • | • |
| A00000356 Magnetic stir bar, Ø 8x40 mm | • | • | • |
| A00001061 Magnetic stir bar, Ø 10x60 mm | • | • | • |
| A00000336 Magnetic cross shape stir bar, Ø 10x5 mm | • | • | • |
| A00000352 Magnetic cross shape stir bar, Ø 20x8 mm | • | • | • |
| A00000354 Magnetic disc stir bar, Ø 10x6 mm | • | • | • |
| A00000268 Temperature Probe AISI 316 Ti Ø 3 mm | • | | |
| A00000349 Temperature probe glass-coated Ø 5 mm | • | | |
| F208B0063 VTF Digital Thermoregulator | • | • | |
| A00000382 Extension for support rod | • | • | • |
| A00000280 Clamp for the probe | • | • | |
| A00000351 Handle for AluBlock removal | • | • | • |
| A00001069 Support rod | • | • | • |

■ AM4 Digital PRO

Multi-position digital magnetic stirrer to manage in parallel four heating plates independently, up to a volume of 20 L each. SpeedServo™ system ensures constant speed even when the viscosity changes. The 135 mm diameter CerAlTop™ plates, made of aluminum alloy coated with ceramic, ensures an outstanding temperature homogeneity over the entire surface, optimum resistance to chemical agents and easy cleaning.

SAFETY AND DURABILITY

The control panel is protected against possible damage caused by liquid spills.

PRECISE THERMOREGULATION

Excellent thermoregulation performance with Pt100 Probe and VTF.



Volume (H₂O)
20 L (x4)



Stirring speed
From 50 to 1500 rpm



Temperature range
From room temperature
to 370 °C




External connection
Pt 100 Probe



External connection
VTF

CONFIGURATIONS - AM4 Digital PRO

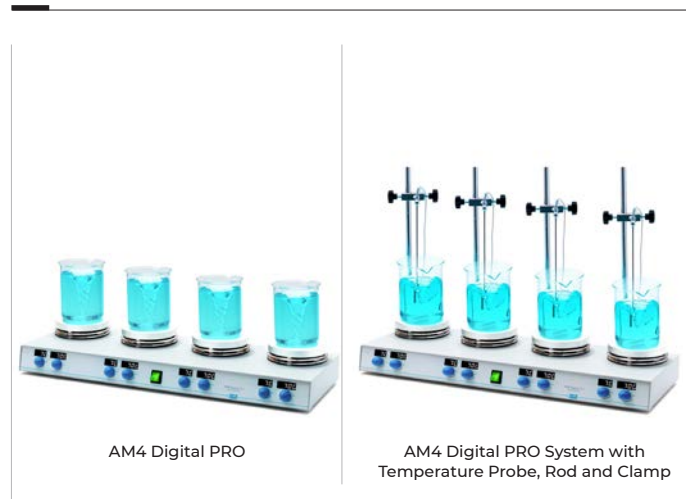

AM4 Digital PRO

AM4 Digital PRO System with
Temperature Probe, Rod and Clamp

INSTRUMENT - CODE

| | | |
|---|------------------|------------|
| AM4 Digital PRO | 230 V / 50-60 Hz | F20500425 |
| AM4 Digital PRO with Temperature Probe, Rod and Clamp | 230 V / 50-60 Hz | SA20500425 |



AM4

Multi-position hot plate stirrer with four separately controlled stirring plates (186 mm between them) that guarantees excellent speed control and resistance to chemicals. SpeedServo™ technology ensures constant speed even when the viscosity changes.

RELIABLE

Aluminum alloy top plate (135 mm diameter) ensures a great temperature homogeneity and optimum heat transfer across the entire surface.



Volume (H₂O)
15 L (x4)



Stirring speed
Up to 1500 rpm



Temperature range
From room temperature
to 370 °C

INSTRUMENT - CODE

AM4

230 V / 50-60 Hz

F20500420



ACCESSORIES - AM4 Series

| | AM4 Digital PRO | AM4 |
|---|-----------------|-----|
| A00001062 Magnetic stir bar, Ø 3x6 mm | • | • |
| A00001063 Magnetic stir bar, Ø 4,5x12 mm | • | • |
| A00001057 Magnetic stir bar, Ø 6x20 mm | • | • |
| A00001056 Magnetic stir bar, Ø 6x35 mm | • | • |
| A00000356 Magnetic stir bar, Ø 8x40 mm | • | • |
| A00001061 Magnetic stir bar, Ø 10x60 mm | • | • |
| A00000336 Magnetic cross shape stir bar, Ø 10x5 mm | • | • |
| A00000352 Magnetic cross shape stir bar, Ø 20x8 mm | • | • |
| A00000354 Magnetic disc stir bar, Ø 10x6 mm | • | • |
| A00000355 Magnetic disc stir bar, Ø 20x10 mm | • | • |
| A00000268 Temperature Probe AISI 316 Ti Ø 3 mm | • | • |
| A00000349 Temperature probe glass-coated Ø 5 mm | • | • |
| F208B0063 VTF Digital Thermoregulator | • | • |
| A00000382 Extension for support rod | • | • |
| A00000280 Clamp for the probe | • | • |
| A00000351 Handle for AluBlock removal | • | • |
| A00001069 Support rod | • | • |

PTFE SAFETY COVER - ACCESSORIES

PTFE safety cover are designed to increase protection and reduce heat loss when working with round bottom flasks and VELP hemispheric bowls.

| | | | |
|---|------------|---|------------|
| PTFE Safety Cover Hemispheric Bowl 50 ml | A00000342* | PTFE Safety Cover Hemispheric Bowl 500 ml | A00000345* |
| PTFE Safety Cover Hemispheric Bowl 100 ml | A00000343* | PTFE Safety Cover Hemispheric Bowl 100 ml | A00000346* |
| PTFE Safety Cover Hemispheric Bowl 250 ml | A00000344* | | |

* not available in USA, Mexico and Canada

HEMISPHERIC BOWLS™ - ACCESSORIES

The Hemispheric Bowls are the most suitable solution to replace oil baths and mantles in order to perform clean and safe analysis with round bottom flasks. VELP Hemispheric Bowl precisely fit the 135 mm diameter plate and the flask surface, ensuring the fastest heat transfer

| | | VOLUME | | | VOLUME |
|---|-----------|--------|---|-----------|--------|
|  | A00000373 | 25 ml |  | A00000333 | 500 ml |
|  | A00000330 | 50 ml |  | A00000334 | 1 L |
|  | A00000331 | 100 ml |  | A00000374 | 3 L |
|  | A00000332 | 250 ml |  | A00000375 | 5 L |



ALUBLOCK™ - ACCESSORIES

The new MonoAluBlock™ and MultiAluBlock™ permits to perform numerous reactions, reducing time and bench space also achieving the highest stirring and heating performances.

MonoAluBlock™ are the ideal solution for all the application that require to stir multiple vials of the same dimensions. MultiAluBlock™ consist of small combinable segments in order to accommodate different sizes of vials at the same time. The user can select the most suitable mix for his application

| ■ MonoAluBlock™ | | NUMBER OF PLACES | VOLUME (ml) | VIAL DIAMETER (mm) | HOLE DEPTH (mm) |
|---|-----------|------------------|-------------|--------------------|-----------------|
|  | A00000340 | 17 | 40 | 28 | 43 |
|  | A00000339 | 17 | 30 | 28 | 30 |
|  | A00000338 | 17 | 20 | 28 | 24 |
|  | A00000341 | 40 | 2 | 12 | 14 |
| ■ MultiAluBlock™ | | NUMBER OF PLACES | VOLUME (ml) | VIAL DIAMETER (mm) | HOLE DEPTH (mm) |
|  | A00000324 | 4 | 40 | 28 | 43 |
|  | A00000325 | 4 | 30 | 28 | 30 |
|  | A00000326 | 4 | 20 | 28 | 24 |
|  | A00000327 | 4 | 16 | 21 | 31 |
|  | A00000328 | 8 | 8 | 17 | 26 |
|  | A00000329 | 11 | 4 | 15 | 20 |
|  | A00000337 | 11 | 2 | 12 | 14 |
|  | A00000323 | 1 | - | - | - |

■ AREC Connect

Innovative, safe and accurate Ceramic Hot Plate Stirrer with embedded Wi-Fi technology to monitor and control laboratory reactions through VELP Ermes cloud platform. Precise thermoregulation and powerful stirring for microbiological, pharmaceutical, chemical and cosmetic laboratories requiring accuracy and efficiency.

OUTSTANDING PERFORMANCE

The brushless motors ensure longevity and guarantee powerful stirring. The cutting-edge heating plate technology ensures temperature homogeneity and stability up to 550 °C. The stirring volume of up to 20 L (H₂O) with excellent coupling.

PRECISE THERMOREGULATION

Excellent thermoregulation performance with Pt100 probe. Thermoregulation range from room temperature to 300 °C with 0.5 °C accuracy.

SAFE AND RESISTANT

The high-quality ceramic top and technopolymer structure ensure premium resistance to chemicals. Safe hot plate warnings when the plate temperature is above 50 °C. The Lock button to prevent unintentional changes

INTUITIVE AND EASY-TO-USE

Easy to operate digital interface with icons.



ermes enabled



Volume (H₂O)
20 L



Stirring speed
From 30 to 1700 rpm




Temperature range
From room temperature
to 550 °C



External connection
Pt 100 Probe

CONFIGURATIONS - AREC Connect

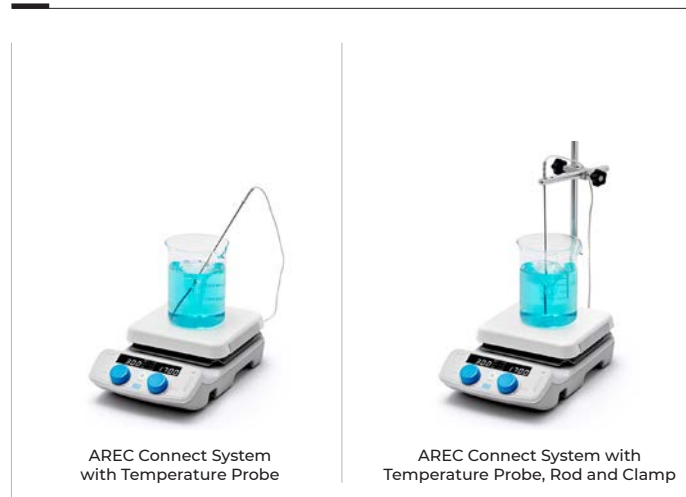

AREC Connect System
with Temperature Probe

AREC Connect System with
Temperature Probe, Rod and Clamp

INSTRUMENT - CODE

| | | |
|--|------------------|------------|
| AREC Connect with Temperature Probe | 230 V / 50-60 Hz | SA20500550 |
| AREC Connect with Temperature Probe | 115 V / 60 Hz | SA20510550 |
| AREC Connect with Temperature Probe, Rod and Clamp | 230 V / 50-60 Hz | SB20500550 |
| AREC Connect with Temperature Probe, Rod and Clamp | 115 V / 60 Hz | SB20510550 |

Included: Magnetic stir bar Ø 8x40 mm.



■ AREC.X

The AREC.X is designed for accurate thermoregulation offering state-of-the-art stirring and heating performance, superior safety features and ease-of-use. The ideal solution for laboratories requiring precise control of the medium's temperature.

EXCELLENT PERFORMANCE

The cutting-edge heating plate technology that ensures temperature homogeneity and stability from room temperature up to 550 °C. Stirring volume of up to 15 L (H₂O).

PRECISE THERMOREGULATION

Optimum thermoregulation of the medium with Pt100 probe (± 1 °C) and VTF Digital Thermoregulator (± 0.5 °C). Thermoregulation range from room temperature to 250 °C.

SAFE AND RESISTANT

The high-quality ceramic top and technopolymer structure ensure premium resistance to chemicals. Safe hot plate warnings when the plate temperature is above 50 °C. The Run-off groove is designed to protect the electronics.



Volume (H₂O)
15 L



Stirring speed
Up to 1500 rpm



Temperature range
From room temperature
to 550 °C




External connection
Pt 100 Probe



External connection
VTF



CONFIGURATIONS - AREC.X

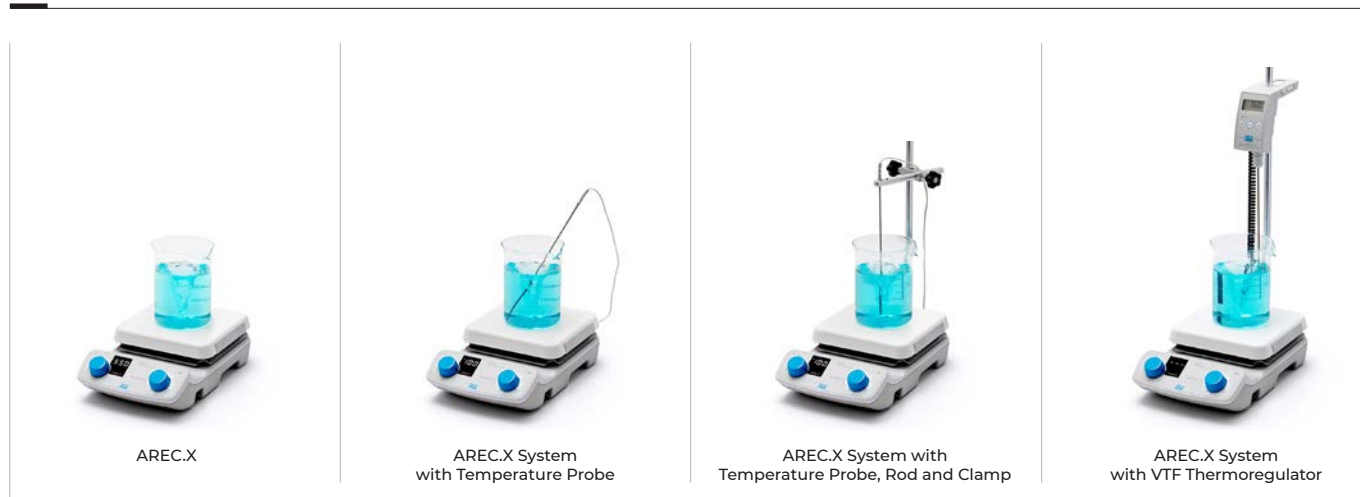

INSTRUMENT - CODE

| | | |
|--|------------------|------------|
| AREC.X | 230 V / 50-60 Hz | F20500554 |
| AREC.X | 115 V / 60 Hz | F20510554 |
| AREC.X with Temperature Probe | 230 V / 50-60 Hz | SA20500554 |
| AREC.X with Temperature Probe | 115 V / 60 Hz | SA20510554 |
| AREC.X with Temperature Probe, Rod and Clamp | 230 V / 50-60 Hz | SB20500554 |
| AREC.X with Temperature Probe, Rod and Clamp | 115 V / 60 Hz | SB20510554 |
| AREC.X with VTF | 230 V / 50-60 Hz | SC20500554 |
| AREC.X with VTF | 115 V / 60 Hz | SC20510554 |

Included: Magnetic stir bar \varnothing 8x40 mm.



■ AREC

The AREC is the safe and reliable hot plate stirrer with ceramic top guaranteeing excellent heat transfer and stirring performance for basic laboratory applications not requiring thermoregulation of the medium.

EXCELLENT PERFORMANCE

The cutting-edge heating plate technology that ensures temperature homogeneity and stability from room temperature up to 550 °C. Stirring volume of up to 15 L (H₂O).

SAFE AND RESISTANT

The high-quality ceramic top and technopolymer structure ensure premium resistance to chemicals. Safe hot plate warnings when the plate temperature is above 50°C. The Run-off groove is designed to protect the electronics.



Volume (H₂O)
15 L



Stirring speed
Up to 1500 rpm



Temperature range
From room temperature
to 550 °C



INSTRUMENT - CODE

| | | |
|------|------------------|-----------|
| AREC | 230 V / 50-60 Hz | F20500551 |
| AREC | 115 V / 60 Hz | F20510551 |

Included: Magnetic stir bar Ø 8x40 mm.



■ HSC

The HSC is the safe and powerful analog hot plate stirrer for applications that do not require accurate temperature control of the medium. The ideal solution for academia's chemistry and biology laboratories.

ENTRY LEVEL SOLUTION

The cutting-edge heating plate technology that ensures temperature homogeneity and stability from room temperature up to 400 °C. Stirring volume of up to 15 L (H₂O).

SAFE AND RESISTANT

The high-quality ceramic top and technopolymer structure ensure premium resistance to chemicals. The Run-off groove is designed to protect the electronics.



Volume (H₂O)
15 L



Stirring speed
Up to 1300 rpm



Temperature range
From room temperature
to 400 °C



INSTRUMENT - CODE

| | | |
|-----|------------------|-----------|
| HSC | 230 V / 50-60 Hz | F20500552 |
| HSC | 115 V / 60 Hz | F20510552 |



ACCESSORIES - AREC & HSC

| | AREC Connect | AREC.X | AREC | HSC |
|--|--------------|--------|------|-----|
| A00000273 Hemispheric bowl for 50 ml flasks | • | • | • | • |
| A00000258 Hemispheric bowl for 100 ml flasks | • | • | • | • |
| A00001071 Hemispheric bowl for 250 ml flasks | • | • | • | • |
| A00001072 Hemispheric bowl for 500 ml flasks | • | • | • | • |
| A00001073 Hemispheric bowl for 1000 ml flasks | • | • | • | • |
| A00001056 Magnetic stir bar, Ø 6x35 mm | • | • | • | • |
| A00000356 Magnetic stir bar, Ø 8x40 mm | • | • | • | • |
| A00000387 Magnetic stir bar, Ø 8x50 mm | • | • | • | • |
| A00001061 Magnetic stir bar, Ø 10x60 mm | • | • | • | • |
| A00000355 Magnetic disc stir bar Ø 20x10mm | • | • | • | • |
| A00001069 Support rod | • | • | • | • |
| A00000382 Extension for support rod | • | • | • | • |
| A00000280 Clamp for PT100 probe | • | • | • | • |
| A00000268 Temperature Probe AISI 316 Ti Ø 3 mm | • | • | • | • |
| A00000349 Temperature Probe glass coated Ø 5 mm | • | • | • | • |
| F208B0063 VTF Digital Thermoregulator | | • | | |
| A00000421 IQ/OQ AREC Connect Manual | • | | | |
| E00010012 VELP Ermes 1 Year Connection | • | | | |
| E00010036 VELP Ermes 3 Years Connection | • | | | |



■ MULTI-HS 15 Digital

Digital multi-position magnetic stirrer with heating designed to guarantee excellent stirring and heating performance for application with 15 positions working under identical conditions. The MULTI-HS presents exclusive features such as timer, autoreverse, safety lock, SpeedServo™ technology and CerAlTop™ plate.

EXCELLENT PERFORMANCE

The brushless motors ensure longevity and guarantee constant speed in all positions.

SAFE AND PROGRAMMABLE

Safe hot plate warnings when the plate temperature is above 50 °C.

Lock feature to prevent unintentional modifications.

Unattended stirring from 1 minute to 99 h : 59 minutes.

EASY TO USE

Easy to operate digital interface with icons. An acoustic signal/sound informs you when the countdown is over.



Volume (H₂O)
250 ml (x15)



Stirring speed
From 100 to 1000 rpm



Temperature range
From room temperature
to 120 °C



Timer &
Auto-reverse



Safety Lock



INSTRUMENT - CODE

MULTI-HS 15 Digital

230 V / 50-60 Hz

F20500531

MULTI-HS 15 Digital

115 V / 60 Hz

F20510531



MULTI-HS 6 Digital

Digital multi-position magnetic stirrer with heating designed to guarantee excellent stirring and heating performance for application with 6 positions working under identical conditions. The MULTI-HS presents exclusive features such as timer, autoreverse, safety lock, SpeedServo™ technology and CerAlTop™ plate.

EXCELLENT PERFORMANCE

The brushless motors ensure longevity and guarantee constant speed in all positions.

SAFE AND PROGRAMMABLE

Safe hot plate warnings when the plate temperature is above 50 °C.

Lock feature to prevent unintentional modifications.

Unattended stirring from 1 minute to 99 h : 59 minutes.

EASY TO USE

Easy to operate digital interface with icons. An acoustic signal/sound informs you when the countdown is over.



Volume (H₂O)
400 ml (x6)



Stirring speed
From 100 to 1000 rpm



Temperature range
From room temperature
to 120 °C



Timer &
Auto-reverse



Safety Lock



INSTRUMENT - CODE

| | | |
|--------------------|------------------|-----------|
| MULTI-HS 6 Digital | 230 V / 50-60 Hz | F20500530 |
| MULTI-HS 6 Digital | 115 V / 60 Hz | F20510530 |



ACCESSORIES - MULTI-HS Series

| | MULTI-HS 6 Digital | MULTI-HS 15 Digital |
|---|--------------------|---------------------|
| A00001056 Magnetic stir bar Ø6x35 mm | • | • |
| A00001057 Magnetic stir bar Ø6x20 mm | • | • |



VTF

The VTF thermoregulator guarantees premium performance in terms of temperature control with an accuracy of ± 0.5 °C. It automatically adjusts the thermoregulation to different factors such as power, load and thermal dispersion.

EASY TO USE

The VTF is equipped with a unique and innovative integrated system to simplify installation on the support rod.



Temperature range
From -10 °C to 300 °C



Timer
From 0 minute
to 24:59 hours



INSTRUMENT - CODE

VTF

12 V dc

F208B0063



ACCESSORIES - VTF

| | VTF |
|--|-----|
| A0000001 Derivation element PW10 | • |
| A0000002 Probe extension cable, 1 m | • |
| A0000003 Glass probe | • |
| A00000280 Clamp for Pt100 Probe | • |
| A00001069 Support rod | • |

Heating Plates

■ RC

Single-position heating plate with aluminum top and temperature regulation. It is designed for the longest lifespan and guarantees excellent resistance to chemicals. It accepts different combinations of AluBlocks™ that can maximize temperature control and save handling time.

SAFETY AND RESISTANCE

Control panel is protected from liquid spills with run-off groove. IP 42 Certified.

EASY TO USE AND UNIVERSAL

Inclined control panel for easy access and visibility. Choice of mains voltage/frequency model.



Volume (H₂O)
15 L



Temperature range
From room temperature
to 370 °C



INSTRUMENT - CODE

| | | |
|----|------------------|-----------|
| RC | 230 V / 50-60 Hz | F20700431 |
| RC | 115 V / 60 Hz | F20710431 |



■ RC2

Two-position heating plate with aluminum top and temperature regulation. It is designed for the longest lifespan and guarantees excellent resistance to chemicals. It accepts different combinations of AluBlocks™ that can maximize temperature control and save handling time.

SAFETY AND RESISTANCE

Control panel shielded against spillage



Volume (H₂O)
15 L (x2)



Temperature range
From room temperature
to 370 °C



INSTRUMENT - CODE

| | | |
|-----|------------------|-----------|
| RC2 | 230 V / 50-60 Hz | F20700430 |
| RC2 | 115 V / 60 Hz | F20710430 |



■ REC

The VELP REC is the laboratory digital heating plate with ceramic top ensuring high resistance to corrosion, and precise temperature setting.

SAFETY AND RESISTANCE

The high-quality ceramic top and technopolymer structure ensure premium resistance to chemicals. The Ceramic heating plate guarantees optimum heat transfer on the entire surface and resistance to chemicals.

EASY TO USE

The bright digital display guarantees a precise temperature setting and showing actual and setpoint temperature. The Hot Surface Safety message displayed when the temperature is above 50 °C.



Volume (H₂O)
15 L



Temperature range
From room temperature
to 550 °C



INSTRUMENT - CODE

| | | |
|-----|------------------|-----------|
| REC | 230 V / 50-60 Hz | F20700553 |
| REC | 115 V / 60 Hz | F20710553 |



ACCESSORIES - HEATING PLATES

| | RC | RC2 | REC |
|---|----|-----|-----|
| A00000373 Hemispheric bowl for 25 ml flasks | • | • | |
| A00000330 Hemispheric bowl for 50 ml flasks | • | • | |
| A00000273 Hemispheric bowl for 50 ml flasks | | | • |
| A00000331 Hemispheric bowl for 100 ml flasks | • | • | |
| A00000258 Hemispheric bowl for 100 ml flasks | | | • |
| A00000332 Hemispheric bowl for 250 ml flasks | • | • | |
| A00001071 Hemispheric bowl for 250 ml flasks | | | • |
| A00000333 Hemispheric bowl for 500 ml flasks | • | • | |
| A00001072 Hemispheric bowl for 500 ml flasks | | | • |
| A00001073 Hemispheric bowl for 1000 ml flasks | | | • |
| A00000334 Hemispheric bowl for 1000 ml flasks | • | • | |
| A00000374 Hemispheric bowl for 3 L flasks | • | • | |
| A00000375 Hemispheric bowl for 5 L flasks | • | • | |
| A00000323 MultiAluBlock™ Base Ø135 mm | • | • | |
| A00000324 MultiAluBlock™, 4 pos. Ø 28 x h 43 mm | • | • | |
| A00000325 MultiAluBlock™, 4 pos. Ø 28 x h 30 mm | • | • | |
| A00000326 MultiAluBlock™, 4 pos. Ø 28 x h 24 mm | • | • | |
| A00000327 MultiAluBlock™, 4 pos. Ø 21 x h 31 mm | • | • | |
| A00000328 MultiAluBlock™, 8 pos. Ø 17 x h 26 mm | • | • | |
| A00000329 MultiAluBlock™, 11 pos. Ø 15 x h 20 mm | • | • | |
| A00000337 MultiAluBlock™, 11 pos. Ø 12 x h 14 mm | • | • | |
| A00000340 MonoAluBlock, 17 pos. Ø 28 x h 43 mm | • | • | |
| A00000339 MonoAluBlock, 17 pos. Ø 28 x h 30 mm | • | • | |
| A00000338 MonoAluBlock, 17 pos. Ø 28 x h 24 mm | • | • | |
| A00000341 MonoAluBlock, 40 pos. Ø 12 x h 14 mm | • | • | |
| A00000342 PTFE Safety cover for bowl 50 ml | • | • | |
| A00000343 PTFE Safety cover for bowl 100 ml | • | • | |
| A00000344 PTFE Safety cover for bowl 250 ml | • | • | |
| A00000345 PTFE Safety cover for bowl 500 ml | • | • | |
| A00000346 PTFE Safety cover for bowl 1000 ml | • | • | |
| A00001069 Support rod | • | • | • |
| A00000382 Extension for support rod | | | • |
| A00000351 Handle for AluBlock removal | • | • | |

Overhead Stirrers

■ OHS 200 Advance

Powerful Wi-Fi Enabled Laboratory overhead stirrer with a TFT 3.5" Display that shows all setting information at a glance. Wi-Fi Connection to VELP Ermes Platform allows remote monitoring and control of operations in real-time. USB interface for use with optional PC software to data log, monitor, turn off/on. Pt100 probe included allowing real-time temperature monitoring.

CONSTANT TORQUE AT ANY SPEED






The powerful brushless motor and electronic control guarantee 200 Ncm torque (speed gear 1) at any speed. Double speed gear for an excellent performance for high volumes and high viscosities.

EASY TO USE

Smartchuck™ system simplifies the shaft accommodation without tools and with one hand only.

SAFE AND RESISTANT

The resistant sealed housing guarantees durability in the most aggressive environments. IP 54 certified.

| | | | | | |
|---|---|---|---|---|---|
|  |  |  |  |  |  |
| Volume (H ₂ O) 100 L | Viscosity Up to 100.000 mPa*s | Torque GEAR 1: 200 Ncm | Stirring speed GEAR 1: 6-400 rpm GEAR 2: 30-2000 rpm | Timer and programmable methods | Safety Lock |



ermes enabled

INSTRUMENT - CODE

| | | |
|-----------------|------------------|-----------|
| OHS 200 Advance | 230 V / 50-60 Hz | F20100480 |
| OHS 200 Advance | 115 V / 60 Hz | F20110480 |



■ OHS 100 Advance

Powerful Wi-Fi Enabled Laboratory overhead stirrer with a TFT 3.5" Display that shows all setting information at a glance. Wi-Fi Connection to VELP Ermes Platform allows remote monitoring and control of operations in real-time. USB interface for use with optional PC software to data log, monitor, turn off/on. Pt100 probe included allowing real-time temperature monitoring.

CONSTANT TORQUE AT ANY SPEED







The powerful brushless motor and electronic control guarantee 100 Ncm torque at any speed.

EASY TO USE

Smartchuck™ system simplifies the shaft accommodation without tools and with one hand only.

SAFE AND RESISTANT

The resistant sealed housing guarantees durability in the most aggressive environments. IP 54 certified.

| | | | | | |
|---|---|---|---|---|---|
|  |  |  |  |  |  |
| Volume (H ₂ O) 100 L | Viscosity Up to 70.000 mPa*s | Torque 100 Ncm | Stirring speed From 30 to 1300 rpm | Timer and programmable methods | Safety Lock |



ermes enabled

INSTRUMENT - CODE

| | | |
|-----------------|------------------|-----------|
| OHS 100 Advance | 230 V / 50-60 Hz | F20100481 |
| OHS 100 Advance | 115 V / 60 Hz | F20110481 |



■ OHS 60 Advance

Powerful Wi-Fi Enabled Laboratory overhead stirrer with a TFT 3.5" Display that shows all setting information at a glance. Wi-Fi Connection to VELP Ermes Platform allows remote monitoring and control of operations in real-time. USB interface for use with optional PC software to data log, monitor, turn off/on. Pt100 probe included allowing real-time temperature monitoring.

CONSTANT TORQUE AT ANY SPEED

The powerful brushless motor and electronic control guarantee 60 Ncm torque at any speed.

EASY TO USE






Smartchuck™ system simplifies the shaft accommodation without tools and with one hand only.

SAFE AND RESISTANT

The resistant sealed housing guarantees durability in the most aggressive environments. IP 54 certified.



ermes enabled

| | | | | | |
|---|---|---|---|---|---|
|  |  |  |  |  |  |
| Volume (H ₂ O) 40 L | Viscosity Up to 50.000 mPa*s | Torque 60 Ncm | Stirring speed From 30 to 2000 rpm | Timer and programmable methods | Safety Lock |

INSTRUMENT - CODE

| | | |
|----------------|------------------|-----------|
| OHS 60 Advance | 230 V / 50-60 Hz | F20100482 |
| OHS 60 Advance | 115 V / 60 Hz | F20110482 |



■ OHS 200 Digital

Powerful overhead stirrer with digital display to precisely set stirring speed and timer. Dedicated 4 segment LED display showing the torque trend. Intuitive knob to navigate the menu and start stop operations.

CONSTANT TORQUE AT ANY SPEED

The powerful brushless motor and electronic control guarantee 200 Ncm torque (speed gear 1) at any speed. Double speed gear for an excellent performance for high volumes and high viscosities.







EASY TO USE

Smartchuck™ system simplifies the shaft accommodation without tools and with one hand only.

SAFE AND RESISTANT

The resistant sealed housing guarantees durability in the most aggressive environments. IP 54 certified.



| | | | | | |
|---|---|---|---|---|---|
|  |  |  |  |  |  |
| Volume (H ₂ O) 100 L | Viscosity Up to 100.000 mPa*s | Torque GEAR 1: 200 Ncm | Stirring speed GEAR 1: 6-400 rpm GEAR 2: 30-2000 rpm | Timer | Safety Lock |

INSTRUMENT - CODE

| | | |
|-----------------|------------------|-----------|
| OHS 200 Digital | 230 V / 50-60 Hz | F20100490 |
| OHS 200 Digital | 115 V / 60 Hz | F20110490 |



■ OHS 100 Digital

Powerful overhead stirrer with digital display to precisely set stirring speed and timer. Dedicated 4 segment LED display showing the torque trend. Intuitive knob to navigate the menu and start stop operations.

CONSTANT TORQUE AT ANY SPEED

The powerful brushless motor and electronic control guarantee 100 Ncm torque at any speed.


EASY TO USE

Smartchuck™ system simplifies the shaft accommodation without tools and with one hand only.

SAFE AND RESISTANT

The resistant sealed housing guarantees durability in the most aggressive environments. IP 54 certified.



| | | | | | |
|---|---|---|---|---|---|
|  |  |  |  |  |  |
| Volume (H ₂ O) 100 L | Viscosity Up to 70.000 mPa*s | Torque 100 Ncm | Stirring speed From 30 to 1300 rpm | Timer | Safety Lock |

INSTRUMENT - CODE

| | | |
|-----------------|------------------|-----------|
| OHS 100 Digital | 230 V / 50-60 Hz | F20100491 |
| OHS 100 Digital | 115 V / 60 Hz | F20110491 |

■ OHS 60 Digital

Powerful overhead stirrer with digital display to precisely set stirring speed and timer. Dedicated 4 segment LED display showing the torque trend. Intuitive knob to navigate the menu and start stop operations.

CONSTANT TORQUE AT ANY SPEED

The powerful brushless motor and electronic control guarantee 60 Ncm torque at any speed.







EASY TO USE

Smartchuck™ system simplifies the shaft accommodation without tools and with one hand only.

SAFE AND RESISTANT

The resistant sealed housing guarantees durability in the most aggressive environments. IP 54 certified.



| | | | | | |
|---|---|---|---|---|---|
|  |  |  |  |  |  |
| Volume (H ₂ O) 40 L | Viscosity Up to 50.000 mPa*s | Torque 60 Ncm | Stirring speed From 30 to 2000 rpm | Timer | Safety Lock |

INSTRUMENT - CODE

| | | |
|----------------|------------------|-----------|
| OHS 60 Digital | 230 V / 50-60 Hz | F20100492 |
| OHS 60 Digital | 115 V / 60 Hz | F20110492 |

■ OHS 40 Digital

Powerful overhead stirrer with digital display to precisely set stirring speed and timer. Dedicated 4 segment LED display showing the torque trend. Intuitive knob to navigate the menu and start stop operations.

CONSTANT TORQUE AT ANY SPEED

The powerful brushless motor and electronic control guarantee 40 Ncm torque at any speed.

EASY TO USE

Smartchuck™ system simplifies the shaft accommodation without tools and with one hand only.

SAFE AND RESISTANT

The resistant sealed housing guarantees durability in the most aggressive environments. IP 54 certified.



Volume (H₂O)
25 L



Viscosity
Up to 25.000 mPa*s



Torque
40 Ncm



Stirring speed
From 30 to 2000 rpm



Timer



Safety Lock

INSTRUMENT - CODE

| | | |
|----------------|------------------|-----------|
| OHS 40 Digital | 230 V / 50-60 Hz | F20100493 |
| OHS 40 Digital | 115 V / 60 Hz | F20110493 |



■ OHS 20 Digital

Powerful overhead stirrer with digital display to precisely set stirring speed and timer. Dedicated 4 segment LED display showing the torque trend. Intuitive knob to navigate the menu and start stop operations.

CONSTANT TORQUE AT ANY SPEED

The powerful brushless motor and electronic control guarantee 20 Ncm torque at any speed.

EASY TO USE

Smartchuck™ system simplifies the shaft accommodation without tools and with one hand only.

SAFE AND RESISTANT

The resistant sealed housing guarantees durability in the most aggressive environments. IP 54 certified.



Volume (H₂O)
25 L



Viscosity
Up to 10.000 mPa*s



Torque
20 Ncm



Stirring speed
From 30 to 2000 rpm



Timer



Safety Lock

INSTRUMENT - CODE

| | | |
|----------------|------------------|-----------|
| OHS 20 Digital | 230 V / 50-60 Hz | F20100494 |
| OHS 20 Digital | 115 V / 60 Hz | F20110494 |



DLH

Digital overhead stirrer suitable for all medium to high-viscosity mixing tasks. The bright display provides easy-to-read details of four important parameters: set speed, actual speed, torque and time.

SAFETY AND EASY TO USE

The DLH is an intuitive overhead stirrer ensuring safety with overload, over-current and overheating protection.

COMPACT AND LIGHTWEIGHT

The special technopolymer structure and the IP 40 Protection Rating guarantee chemical resistance and easy handling, being among the lightest overhead stirrer on the market.



Volume (H₂O)
40 L



Viscosity
Up to 50.000 mPa*s



Torque
80 Ncm



Stirring speed
From 50 to 2000 rpm



Timer

INSTRUMENT - CODE

DLH

110-230 V / 50-60 Hz

F201A0157



DLS

Digital overhead stirrer suitable for all medium viscosity mixing tasks. The bright display provides easy-to-read details of four important parameters: set speed, actual speed, torque and time.

SAFETY AND EASY TO USE

The DLS is an intuitive overhead stirrer ensuring safety with overload, over-current and overheating protection.

COMPACT AND LIGHTWEIGHT

The special technopolymer structure and the IP 40 Protection Rating guarantee chemical resistance and easy handling, being among the lightest overhead stirrer on the market.



Volume (H₂O)
25 L



Viscosity
Up to 25.000 mPa*s



Torque
40 Ncm



Stirring speed
From 50 to 2000 rpm



Timer

INSTRUMENT - CODE

DLS

110-230 V / 50-60 Hz

F201A0155

DLS - System

110-230 V / 50-60 Hz

SA201A0155*

* Universal H-Stand, Double Clamp, Stirring Shaft with Propeller included.



■ PW

Laboratory overhead stirrer with analog control suitable for high viscosity mixing applications. Easy to install and to operate thanks to the hand-tightened chuck requiring no tools.

SAFETY AND EASY TO USE

The PW ensures safety with overload, over-current and overheating protection.

COMPACT AND LIGHTWEIGHT

The special technopolymer structure and the IP 40 Protection Rating guarantee chemical resistance and easy handling, being among the lightest overhead stirrer on the market.



Volume (H₂O)
70 L



Viscosity
Up to 100.000 mPa*s



Torque
120 Ncm



Stirring speed
From 20 to 1200 rpm

INSTRUMENT - CODE

PW

110-230 V / 50-60 Hz

F201A0150



■ LH

Laboratory overhead stirrer with analog control suitable for high viscosity mixing applications. Easy to install and to operate thanks to the hand-tightened chuck requiring no tools.

SAFETY AND EASY TO USE

The LH ensures safety with overload, over-current and overheating protection.

COMPACT AND LIGHTWEIGHT

The special technopolymer structure and the IP 40 Protection Rating guarantee chemical resistance and easy handling, being among the lightest overhead stirrer on the market.



Volume (H₂O)
40 L



Viscosity
Up to 50.000 mPa*s



Torque
80 Ncm



Stirring speed
From 50 to 2000 rpm

INSTRUMENT - CODE

LH

110-230 V / 50-60 Hz

F201A0156



■ LS

Laboratory overhead stirrer with analog control suitable for high viscosity mixing applications. Easy to install and to operate thanks to the hand-tightened chuck requiring no tools.

SAFETY AND EASY TO USE

The LS ensures safety with overload, over-current and overheating protection.

COMPACT AND LIGHTWEIGHT

The special technopolymer structure and the IP 40 Protection Rating guarantee chemical resistance and easy handling, being among the lightest overhead stirrer on the market.



Volume (H₂O)
25 L



Viscosity
Up to 25.000 mPa*s



Torque
40 Ncm



Stirring speed
From 50 to 2000 rpm

INSTRUMENT - CODE

| | | |
|-------------|----------------------|-------------|
| LS | 110-230 V / 50-60 Hz | F201A0151 |
| LS - System | 110-230 V / 50-60 Hz | SA201A0151* |

* Universal H-Stand, Double Clamp, Stirring Shaft with Fixed Blade included.



■ ES

Laboratory overhead stirrer with analog control suitable for low viscosity mixing applications. Easy to install and to operate thanks to the hand-tightened chuck requiring no tools.

SAFETY AND EASY TO USE

The ES ensures safety with overload, over-current and overheating protection.

COMPACT AND LIGHTWEIGHT

The special technopolymer structure and the IP 40 Protection Rating guarantee chemical resistance and easy handling, being among the lightest overhead stirrer on the market.



Volume (H₂O)
15 L



Viscosity
Up to 1.000 mPa*s



Torque
15 Ncm



Stirring speed
From 50 to 1300 rpm

INSTRUMENT - CODE

| | | |
|----|----------------------|-----------|
| ES | 110-230 V / 50-60 Hz | F201A0152 |
|----|----------------------|-----------|



ACCESSORIES - OVERHEAD STIRRERS

| | OHS Advance | OHS Digital | DLH | DLS | PW | LH | LS | ES |
|--|-------------|-------------|-----|-----|----|----|----|----|
| A00001301 Double clamp | • | • | • | • | • | • | • | • |
| A00001302 Ribbon clamp | • | • | • | • | • | • | • | • |
| A00001312 Chuck guard OHS | • | • | | | | | | |
| A00000372 Support rod for probe | • | | | | | | | |
| A00000002 Probe extension cable, length 1 m | • | | | | | | | |
| A00000003 Glass probe for VTF | • | | | | | | | |
| TA00000394 Kit Pt100 probe XL for OHS | • | | | | | | | |
| A00000381 IQ/OQ Manual OHS Digital | | • | | | | | | |
| A00000380 IQ/OQ Manual OHS Advance | • | | | | | | | |
| A00000391 OHS/ControllerSoft | • | | | | | | | |
| E00010012 VELP Ermes 1 Year Connection | • | | | | | | | |
| E00010036 VELP Ermes 3 Year Connection | • | | | | | | | |

STIRRING SHAFTS - ACCESSORIES

VELP H-Stands are the ideal solution when the overhead stirrer is working on laboratory benches.

We have designed our H-Stand to guarantee the maximum stability thanks to the H-shape of the base and the maximum resistance choosing premium materials such as aluminum and stainless steel, properly treated in order to resist to chemicals and corrosion.

The Telescopic H-Stand facilitates the daily operations being easy-to-use when raising and lowering the overhead stirrer from your sample.

UNIVERSAL H-STAND



Code A00000366

TELESCOPIC H-STAND



Code A00000369


STIRRING SHAFTS - ACCESSORIES

Overhead stirrers require stirring shaft for precise working. Choosing the correct stirring shaft is an important step to obtain the best stirring performance.


It is important to evaluate the characteristics and the effects on the fluid of the stirring shaft, together with the volume and viscosity:

- Desired type of flow: Axial, Tangential, Radial;
- Speed range of your application.


STIRRING SHAFT WITH FLOATING BLADES (VL - L)

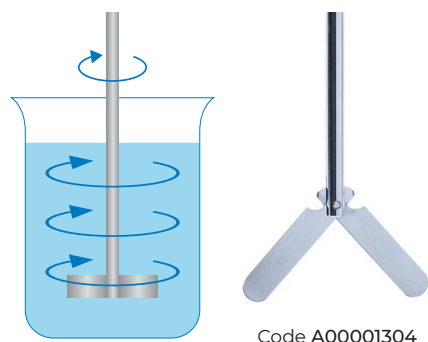


STIRRING SHAFT WITH 6-HOLE PADDLE (L - M)

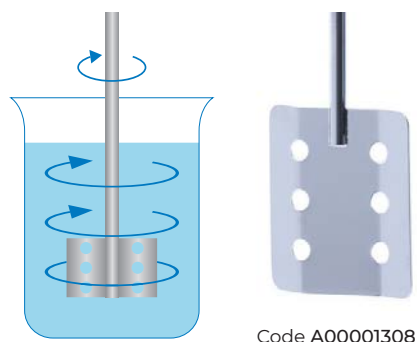


STIRRING SHAFT WITH FIXED BLADE (VL - L - M)

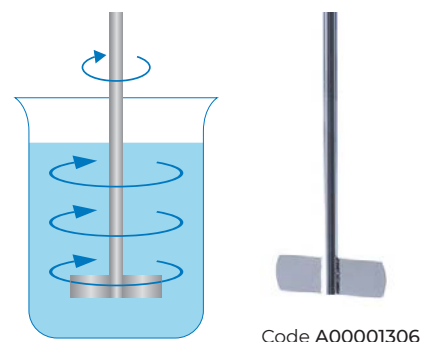


STIRRING SHAFT WITH TURBO PROPELLER (M - H)

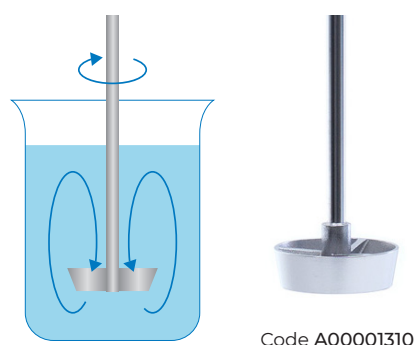


STIRRING SHAFT WITH FOLDING BLADE (VL - L)

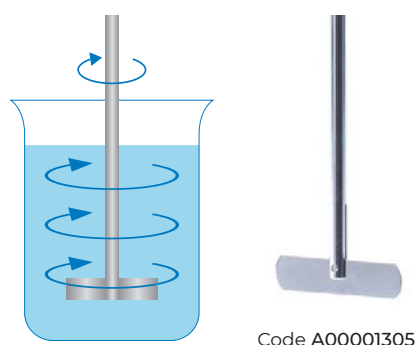


STIRRING SHAFT WITH TURBINE (M - H)

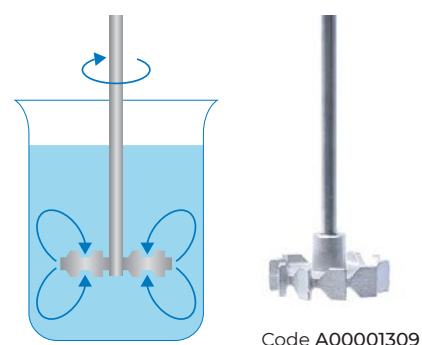

STIRRING SHAFT WITH PROPELLER (VL - L - M)

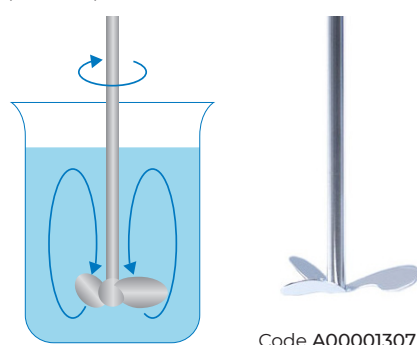

STIRRING SHAFT WITH ANCHOR (M - H)

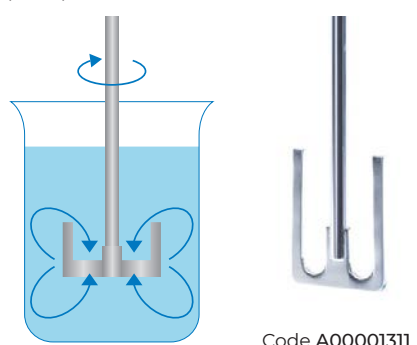

VISCOSITY RANGE

| Viscosity Range | |
|-----------------|------------------|
| Very Low (VL) | 0 - 100 |
| Low (L) | 100 - 1,000 |
| Medium (M) | 1,000 - 10,000 |
| High (H) | 10,000 - 100,000 |

SAMPLE VISCOSITY

| Sample | mPa*s |
|-----------------------|---------|
| Water | 1 |
| Kerosene | 10 |
| Lubricating oil | 100 |
| Castrol oil, Glycerin | 1,000 |
| Shampoo | 3,000 |
| Refined honey | 7,000 |
| Body cream | 8,000 |
| Honey | 10,000 |
| Chocolat syrup | 25,000 |
| Ketchup | 50,000 |
| Toothpaste | 70,000 |
| Molasses, Asphalt | 100,000 |

■ CONTROLLER Advance (not available in USA, Mexico and Canada)

Powerful and innovative solution designed to control brushless motors of Reaction Systems of up to 100 L volume. The Controller Advance features a TFT 3.5" Display that shows all setting information at a glance. Wi-Fi Connection to Velp Ermes Platform to monitor status and data logging on the cloud. USB interface for use with optional PC software to data log, monitor, turn off/on. USB also allow Software Updates.

- Timer for unattended operations
- Real time torque, speed and temperature graph monitoring
- Intuitive knob to navigate the menu and start stop operations
- Programmable methods to setup to 4 Time/Speed Ramps.



INSTRUMENT - CODE

Controller Advance

115-230 V / 50-60 Hz

F201A0500



■ CONTROLLER Digital (not available in USA, Mexico and Canada)

Powerful and innovative solution designed to control brushless motors of Reaction Systems of up to 100 L volume. The Controller Digital features a digital display to set the speed and a 4 segment LED display to monitor torque.

- Timer for unattended operations
- Intuitive knob to navigate the menu and start stop operations
- Clockwise and counterclockwise direction changes.



INSTRUMENT - CODE

Controller Digital

115-230 V / 50-60 Hz

F201A0505



ACCESSORIES - CONTROLLER

| | CONTROLLER ADVANCE | CONTROLLER DIGITAL |
|--|--------------------|--------------------|
| A00000391 OHS/ControllerSoft | • | |
| A00000394 Pt100 Probe XL Ø6 500mm | • | |
| A00000415 Motor 10 lb*in, 471 rpm , 1/11 HP - 24V | • | • |
| A00000416 Motor 28 lb*in, 500 rpm , 1/11 HP - 24V | • | • |
| A00000417 Motor 55 lb*in, 250 rpm , 1/11 HP - 24V | • | • |
| E00010012 VELP Ermes 1 year Connection | • | |
| E00010036 VELP Ermes 3 year Connection | • | |

Vortex Mixers

■ Wizard IR

Vortex mixer featuring the highly innovative and patented IR technology to reduce repetitive strains. IR sensor and Continuous mode available. Large choice of complementary accessories to accommodate tubes, plates and flasks.

EXCLUSIVE ERGONOMICS

The VELP IR sensor uses an infrared system (IR) to detect the presence of the tube and automatically starts vibrating. No pressure is required to maintain vibration, drastically reducing arm stress.

STABLE AND ROBUST

Robust low profile design provides exceptional stability. The technopolymer housing and zinc base ensure optimum chemical resistance and excellent stability on the bench.



Stirring speed
Up to 3000 rpm



Patented
IR Sensor



INSTRUMENT - CODE

Wizard

100÷240 V / 50-60 Hz

F202A0175



■ Classic

Compact vortex mixer suitable for different mixing requirements, with touch or continuous mode. It is highly resistant to chemicals and has 3 anti-slip rubber feet for stability.

ROBUST AND CONVENIENT

The IP 42 Protection Rating and the technopolymer casing guarantee good chemical resistance against particles and liquids. It uses universal voltage/frequency.

VERSATILE

Large choice of accessories to suit tubes, plates and flasks. Accessories fit by simple pressure.



Stirring speed
Up to 3000 rpm



INSTRUMENT - CODE

Classic

100÷240 V / 50-60 Hz

F202A0173



TX4

Exclusive digital vortex mixer featuring the highly innovative and patented IR technology to reduce repetitive strains. The bright digital display ensures precision and reproducibility. The integrated timer allows programming the vortex for unattended operations. Large choice of complementary accessories to accommodate tubes, plates and flasks.

EXCLUSIVE ERGONOMICS

The VELP IR sensor uses an infrared system (IR) to detect the presence of the tube and automatically starts vibrating. No pressure is required to maintain vibration, drastically reducing arm stress.

STABLE AND ROBUST

The technopolymer housing and zinc base ensure optimum chemical resistance and excellent stability on the bench.



Stirring speed
Up to 3000 rpm



Patented
IR Sensor



Timer



INSTRUMENT - CODE

TX4

100÷240 V / 50-60 Hz

F202A0270



ZX4

ZX4 Vortex mixer featuring the highly innovative and patented IR technology to reduce repetitive strains. IR sensor and Continuous mode available. Large choice of complementary accessories to accommodate tubes, plates and flasks.

EXCLUSIVE ERGONOMICS

The VELP IR sensor uses an infrared system (IR) to detect the presence of the tube and automatically starts vibrating. No pressure is required to maintain vibration, drastically reducing arm stress.

STABLE AND ROBUST

The technopolymer housing and zinc base ensure optimum chemical resistance and excellent stability on the bench.



Stirring speed
Up to 3000 rpm



Patented
IR Sensor



INSTRUMENT - CODE

ZX4

100÷240 V / 50-60 Hz

F202A0280

ZX4 and Foam Stands

80÷260 V / 50-60 Hz

SA202A0280*

* Stand for 19 microvials Eppendorf 1,5 ml, Customizable soft foam top and Small rubber supporting plate Ø 50mm included



ZX3

Reliable vortex mixer with adjustable stirring speed, touch or continuous mode. ZX3 features large, anti-slip rubber feet to ensure maximum stability. It is also highly resistant to chemicals.

ROBUST AND CONVENIENT

The IP 42 Protection Rating and the technopolymer casing guarantee good chemical resistance against particles and liquids. It uses universal voltage/frequency.

VERSATILE

Large choice of accessories to suit tubes, plates and flasks. Accessories fit by simple pressure.



Stirring speed
Up to 3000 rpm



INSTRUMENT - CODE

| | | |
|---------------------|----------------------|-------------|
| ZX3 | 100÷240 V / 50-60 Hz | F202A0176 |
| ZX3 and Foam Stands | 100÷240 V / 50-60 Hz | SA202A0176* |

* Stand for 19 microvials Eppendorf 1,5 ml, Customizable soft foam top and Small rubber supporting plate Ø 50mm included



RX3

Basic vortex mixer with fixed stirring speed and touch mode. RX3 features large, anti-slip rubber feet to ensure maximum stability. It is also highly resistant to chemicals.

ROBUST AND CONVENIENT

The IP 42 Protection Rating and the technopolymer casing guarantee high-grade chemical resistance.



Stirring speed
3000 rpm



INSTRUMENT - CODE

| | | |
|-----|----------------------|-----------|
| RX3 | 100÷240 V / 50-60 Hz | F202A0171 |
|-----|----------------------|-----------|



ACCESSORIES - VORTEX MIXERS

| | Wizard | Classic | TX4 | ZX4 | ZX3 | RX3 |
|--|--------|---------|-----|-----|-----|-----|
| A00000012 Foam stand for 19 microvials Eppendorf, 1.5 ml | • | • | • | • | • | |
| A00000013 Customizable soft foam top | • | • | • | • | • | |
| A00000014 Foam stand for 5 test tubes Ø 16 mm | • | • | • | • | • | |
| A00000015 Foam stand for microtiter | • | • | • | • | • | |
| A00000016 Small rubber supporting plate Ø 50 mm | • | • | • | • | • | • |
| A00000017 Rubber supporting plate Ø 94 mm | | • | | | • | |
| A00000019 Foam stand for 4 tubes Ø 29 mm | • | • | • | • | • | |
| A00000389 Foam stand for 5 test-tubes Ø 24 mm | • | • | • | • | • | |
| 10002021 Mixing cup head | • | • | | | | |
| 10005195 Mixing cup head | | | • | • | • | • |
| A000003083 Plug USA for power supply | • | • | • | • | • | • |
| 10003084 Plug UK for power supply | • | • | • | • | • | • |
| 10003085 Plug AU for power supply | • | • | • | • | • | • |

ACCESSORIES - FOAM STAND

VELP vortex mixers can be connected to the dedicated platforms ensuring versatility for a variety of applications. VELP accessories are easy to fix on the vortex with a simple pressure and ensure excellent mixing performance of different types of tubes and containers.

CUSTOMIZABLE SOFT FOAM TOP

Code A00000013



FOAM STAND FOR 5 TEST-TUBES Ø 16 mm

Code A00000014



FOAM STAND FOR 5 TEST-TUBES Ø 24 mm

Code A00000389



FOAM STAND FOR 4 TEST-TUBES Ø 29 mm

Code A00000019



FOAM STAND FOR N° 19 MICROVIALS 1,5 ml

Code A00000012



FOAM STAND FOR MICROTITER

Code A00000015



MULTI-TX5

The MULTI-TX5 Digital is the exclusive solution for all the multi-tube mixing tasks, even the most challenging ones. The MULTI-TX5 Digital supports various foam racks for different vortexing tubes being an extremely versatile multi-tube vortex mixer to emulsify, disintegrate, mix, or homogenize samples.

SAFE AND PROGRAMMABLE

Lock feature to prevent unintentional modifications.
Unattended stirring from 1 minute to 99 h : 59 minutes.

EASY TO USE

Easy to operate digital interface with icons. An acoustic signal/sound informs you when the countdown is over.

SIMPLIFIED OPERATION

With the MULTI-TX5 Digital you can perform any mixing tasks and then change the foam rack for your next application in just a few simple steps.



Stirring speed
From 100 to 2500 rpm



Timer &
Auto-reverse



Safety Lock



INSTRUMENT - CODE

MULTI-TX5 Digital Vortex

100÷240 V / 50-60 Hz

F202A0460



ACCESSORIES - SAMPLE RACK

Various foam racks for vortexing tubes with different dimensions ensure premium versatility in the lab and make the MULTI-TX5 Digital ideal for a wide variety of mixing applications.

SAMPLE RACK FOR TEST TUBE Ø 10 mm, 44 POSITIONS

Code A00000449



SAMPLE RACK FOR TEST TUBE Ø 12 mm, 44 POSITIONS

Code A00000450



SAMPLE RACK FOR TEST TUBE Ø 13 mm, 44 POSITIONS

Code A00000451



SAMPLE RACK FOR TEST TUBE Ø 16 mm, 44 POSITIONS

Code A00000452



SAMPLE RACK FOR TEST TUBE Ø 25 mm, 14 POSITIONS

Code A00000453



SAMPLE RACK FOR TEST TUBE Ø 29 mm, 14 POSITIONS

Code A00000454



Homogenizer

■ OV5

Versatile homogenizer ideal for dispersing, homogenizing and mixing biological tissues, pharmaceutical products, cosmetics and food products. Wide selection of rotors and stators to match the application.

ROBUST AND CONVENIENT

The stainless steel shafts, with PTFE seals, can be easily and quickly connected through a snap hook. OV5 has a smooth start and safety switch.

VERSATILE

OV5 has a wide selection of rotors and stators to match the application.



Volume (H₂O)
8 L



Stirring speed
Up to 3000 rpm

OPTIONAL - ACCESSORIES

| | |
|-----------------------------|-----------|
| Stator rotor release kit | A00000279 |
| Double Rod Mini H-Stand Set | A00000383 |
| Double clamp | A00001301 |
| Strap clamp | A00000044 |



INSTRUMENT - CODE

| | | |
|--------------|---------------|-------------|
| OV5 | 230 V / 50 Hz | R20900010 |
| OV5 - System | 230 V / 50 Hz | SA20900010* |

* Double Rod Mini H-Stand Set, VSS2CSR2 Dispersing tool included



CHOOSE THE MOST SUITABLE DISPERSING TOOL - COMPOSITION

HOMOGENIZER



SHAFT




STATOR



ROTOR



POSSIBLE COMBINATIONS - COMPOSITION

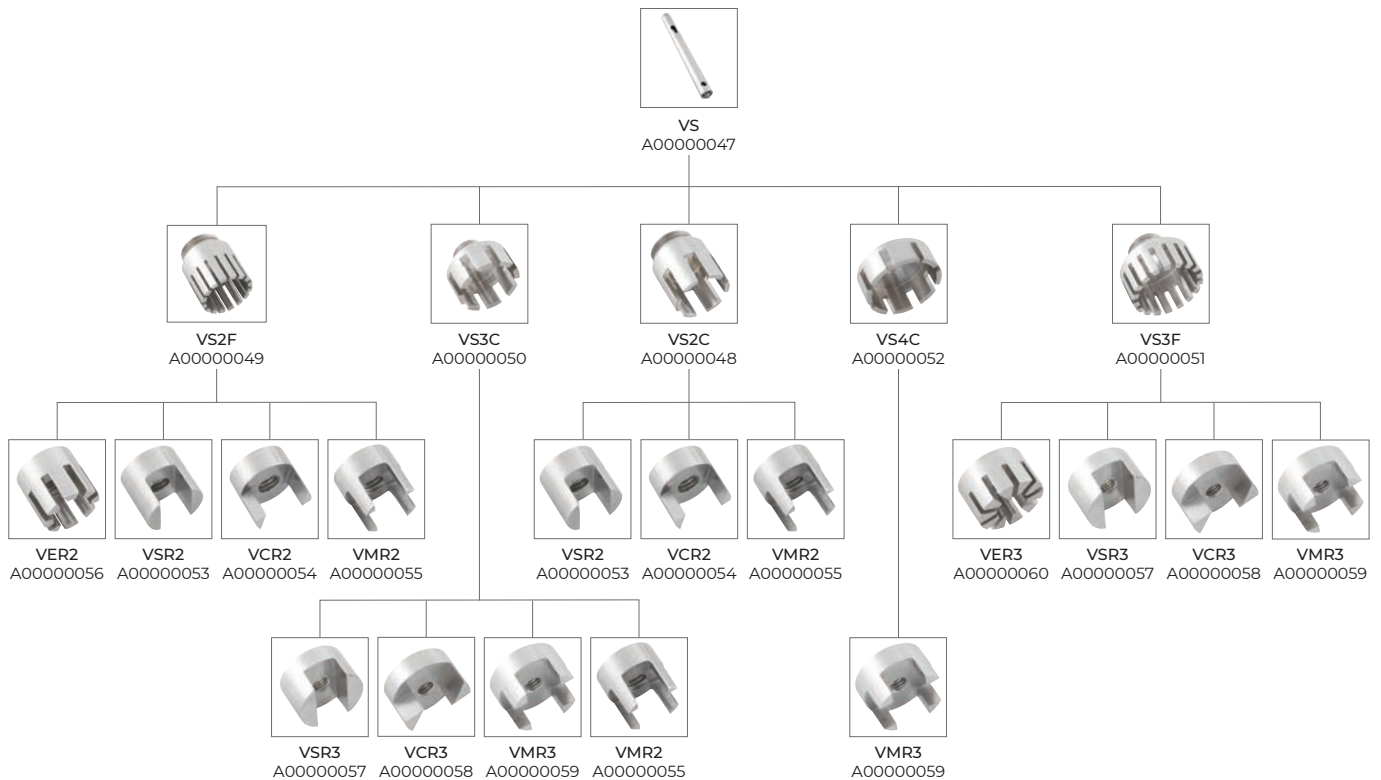

A STAND ALONE SOLUTION FOR A WIDE RANGE OF APPLICATIONS

OV5 Homogenizer is the right solution for dispersing, homogenizing and mixing biological tissue samples (cells, animal and vegetal tissues), pharmaceutical products, cosmetics and food products. It is a versatile instrument: the same model of shaft suitable for all the applications and configurable with different models of rotor and stator, for the specific application of the customer.

POSSIBLE COMBINATIONS - AVAILABLE CONFIGURATIONS

| | Code No | Application field | Function | Treatable volume (water) | Max circum. speed | Ø Rotor | Ø Stator | Tool length | Min/max immersion depth | Ultimate fineness µm Suspension emulsion | |
|------------|-----------|------------------------|----------|--------------------------|-------------------|---------|----------|-------------|-------------------------|--|------|
| * VSS2CSR2 | A00000026 | CE, IF, PC, SI | A | 100-5000 | 22,7 | 15 mm | 20 mm | 220 mm | 40/175 mm | 10-50 | 1-10 |
| VSS2CCR2 | A00000027 | CT, IA, M, SI | B | 100-5000 | 22,7 | 15 mm | 20 mm | 220 mm | 40/175 mm | 10-50 | 1-10 |
| VSS2CMR2 | A00000028 | CE, VE | A | 100-5000 | 22,7 | 15 mm | 20 mm | 220 mm | 40/175 mm | 10-50 | 1-10 |
| *VSS2FER2 | A00000029 | CT, IF, SI, VE | C | 100-5000 | 22,7 | 15 mm | 20 mm | 220 mm | 40/175 mm | 10-50 | 1-10 |
| VSS2FCR2 | A00000031 | BT, CT, IA, IT, M, SI | B | 100-5000 | 22,7 | 15 mm | 20 mm | 220 mm | 40/175 mm | 10-150 | 1-10 |
| VSS2FMR2 | A00000032 | CE, CT, IA, IC, PC, VE | A | 100-5000 | 22,7 | 15 mm | 20 mm | 220 mm | 40/175 mm | 10-150 | 1-10 |
| *VSS3CSR3 | A00000033 | CT, IA, IF, M, SI | A | 100-8000 | 34,9 | 23 mm | 30 mm | 220 mm | 40/175 mm | 5-25 | 1-5 |
| VSS3CCR3 | A00000034 | CT, IA, IF, M, SI | B | 100-8000 | 34,9 | 23 mm | 30 mm | 220 mm | 40/175 mm | 5-25 | 1-5 |
| VSS3CMR3 | A00000035 | CE, VE | A | 100-8000 | 34,9 | 23 mm | 30 mm | 220 mm | 40/175 mm | 5-25 | 1-5 |
| VSS3CMR2 | A00000036 | CE, IA, SI | D | 250-2000 | 34,9 | 15 mm | 30 mm | 220 mm | 40/175 mm | High speed mixer | |
| *VSS3FER3 | A00000037 | CT, IF, SI, VE | C | 100-8000 | 34,9 | 23 mm | 30 mm | 220 mm | 40/175 mm | 5-25 | 1-5 |
| VSS3FSR3 | A00000038 | CT, IF, SI, VE | A | 100-8000 | 34,9 | 23 mm | 30 mm | 220 mm | 40/175 mm | 5-25 | 1-5 |
| VSS3FMR3 | A00000040 | CE, IA, IC, IF, IT | A | 100-8000 | 34,9 | 23 mm | 30 mm | 220 mm | 40/175 mm | 5-25 | 1-5 |
| *VSS4CMR3 | A00000041 | CE, IA, SI | D | 1000-4000 | 34,9 | 23 mm | 40 mm | 220 mm | 40/175 mm | High speed mixer | |
| **VSS5CSR4 | A00000046 | BT, M | A | 0,2-50 | 6,3 | 4 mm | 5 mm | 128 mm | 10/60 mm | 10-50 | 1-10 |

* Most used model

** The dispersing tool works with Ø 4 mm rotor and Ø 5 mm stator for microbiological applications (e.g. suitable for Eppendorf, cuvettes, etc.)

*** BT = biotechnology, CE = ceramic industry, CT = paper & tissue industry, IA = food industry, IC = cosmetic industry, IF = pharmaceutical industry, IT = tobacco industry, M = medicine, PC = petrochemistry industry, SI = sewage pollution control, VE = paint industry

**** A = dispersing tool for solid/liquid media, B = dispersing tool with blades for fibrous/stringy materials, C = dispersing tool for water/oil or oil/water emulsions, D = stirring shaft

COD Thermoreactors

■ ECO 16

Thermoreactor with 16 positions (14 pos. with \varnothing 16 mm and 2 pos. \varnothing 22 mm) that provide easy and intuitive programming of both time and temperature. The unit automatically shuts down when the analysis time expires. The aluminum heating block offers optimum thermal conditions and a high level of homogeneity at all temperatures

INTUITIVE AND DIGITAL

ECO 16 ensures a simple and easy programming of time (from 1 to 199 minutes or continuous) and temperature (from ambient to 160 °C).

FAST AND SAFE

ECO 16 is suitable for performing a short time COD analysis, in only 30 minutes with a high thermoreaction temperature (160 °C). It has also the transparent test tube cover that offers a constant visual control.



Temperature range
From room temperature
to 160 °C



Timer



INSTRUMENT - CODE

| | | |
|--------|------------------|-----------|
| ECO 16 | 230 V / 50-60 Hz | F10100126 |
| ECO 16 | 115 V / 50-60 Hz | F10110126 |



■ ECO 6

Intuitive and digital 6-position thermoreactor for Chemical Oxygen Demand (COD) with settable temperature and timer (up to 199 minutes). The unit automatically shuts down when the analysis time expires. The aluminum heating block offers optimum thermal conditions and a high level of homogeneity at all temperatures

INTUITIVE AND DIGITAL

ECO 6 ensures a simple programming of time (from 1 to 199 minutes or continuous) and temperature (from ambient to 200 °C).

VERSATILE

ECO 6 COD Thermoreactor has been originally designed to support \varnothing 42 mm tubes processing up to 6 x 200 ml, but thanks to smart reducer it can be used with 18 small tubes of 16 mm or 6 x 22 mm too.



Temperature range
From room temperature
to 200 °C



Timer



INSTRUMENT - CODE

| | | |
|-------|------------------|-----------|
| ECO 6 | 230 V / 50-60 Hz | F10100120 |
| ECO 6 | 115 V / 50-60 Hz | F10110120 |



■ ECO 25

Thermoreactor with 25 positions for Ø16mm vials that provide easy and intuitive programming of both time and temperature. The unit automatically shuts down when the analysis time expires. The aluminum heating block offers optimum thermal conditions and a high level of homogeneity at all temperatures

INTUITIVE AND DIGITAL

ECO 25 ensure a simple programming of time, selectable between four different options (30, 60, 120 minutes or continuous) and five temperatures (70, 100, 120, 150 and 160 °C).

SAFE

ECO 25 have a transparent test tube cover and various temperature safety features. The test tube extractor simplifies tubes handling.



Temperature range
70, 100, 120, 150 & 160



Timer



INSTRUMENT - CODE

ECO 25

115-230 V / 50-60 Hz

F101A0125



■ ECO 8

Thermoreactor with 8 positions for Ø16mm vials (plus 1 for Ø22mm vials) providing easy and intuitive programming of both time and temperature. The unit automatically shuts down when the analysis time expires. The aluminum heating block offers optimum thermal conditions and a high level of homogeneity at all temperatures

INTUITIVE AND DIGITAL

ECO 8 ensure a simple programming of time, selectable between four different options (30, 60, 120 minutes or continuous) and five temperatures (70, 100, 120, 150 and 160 °C).

SAFE

ECO 8 have a transparent test tube cover and various temperature safety features.



Temperature range
70, 100, 120, 150 & 160



Timer



INSTRUMENT - CODE

ECO 8

115-230 V / 50-60 Hz

F101A0127



ACCESSORIES - COD THERMOREACTORS

| | ECO 6 | ECO 16 | ECO 8 | ECO 25 |
|---|-------|--------|-------|--------|
| A0000145 COD test tubes Ø 42x200 mm, 200 ml with cone NS 29/32, 3 pcs/box | • | | | |
| A00001041 Air condenser with ground cone | • | | | |
| A00001042 PTFE sheaths for 29/32 cones | • | | | |
| A00001043 Stainless steel handle for moving 6 test tubes Ø 42 mm simultaneously | • | | | |
| A00001045 Antisplash bell | • | | | |
| A00000362 Adapter for probe Ø 3 mm | • | • | • | • |
| A00001044 Anticorrosion reducer Ø 42 mm with 3 holes Ø 16 mm | • | | | |
| A00001046 Anticorrosion reducer Ø 42 mm with 1 hole Ø 22 mm | • | | | |
| A00001051 Safety shield | | • | | |
| A00001039 ECO 25 Test tube extractor | | | | • |
| 10000002 Absorption attachment for condenser NS 29/32 for sludge analysis | | • | • | • |
| CA0091636 Holder for 12 round glass cells | | • | • | • |
| CA0091666 Test tube for sample decomposition Ø 22mm, NS 19/26 cone with glass cap for sludge analysis | | • | • | • |
| CA0091667 Condenser 200 mm type KS with 3-meter polyethylene tube for sludge analysis | | • | • | • |
| CM0091680 Set of 20 glass test tubes with internal Ø 14 mm and external Ø 16 mm | | • | • | • |

BOD and Respirometers

RESPIROMETRIC Sensor System 6 BOD

The RESPIROMETRIC Sensor System BOD is designed to be a versatile solution with minimum footprint for aerobic degradation tests for Industrial and Domestic water quality testing.

- Ideal for: BOD5, BOD7, BOD Ultimate, Assessment of disinfection on BOD.
- Biodegradation of monomers application to study the microbial activities involved in the production of plastics.
- Results reading in Pressure (millibar) and BOD (mg/l).
- Easy-to-use and powerful sensor with 4 digit LED Display.
- Unique wireless control from your PC and VELP Ermes cloud platform.



Wireless DataBox™
and RESPIROSoft™



BOD Scales
90, 250, 600, 999, 4000 ppm

INSTRUMENT - CODE

| | | |
|---|------------------|------------|
| RESPIROMETRIC Sensor System 6 | 230 V / 50-60 Hz | SA10200146 |
| RESPIROMETRIC Sensor System 6 | 115 V / 60 Hz | SA10210146 |
| RESPIROMETRIC Sensor System 6 no Databox™ | 230 V / 50-60 Hz | SB10200146 |
| RESPIROMETRIC Sensor System 6 no Databox™ | 115 V / 60 Hz | SB10210146 |

Included in the system: Stirring station, 500 ml amber bottles, KOH containers, stirring bars, DataBox™ and RESPIROSoft™ software.

ermes enabled



RESPIROMETRIC Sensor System 6 Maxi BMP

The RESPIROMETRIC Sensor System Maxi BMP is designed to support anaerobic respiration studies to determine the Biochemical Methane Potential (BMP) of biomass. The BMP analysis studies the production of biogas from organic wastes mainly for Waste and Wastewater treatment and for Renewable Energy Production.

- Results reading in Pressure (millibar) and automatic conversion in BMP value (NmLCH₄/gSV)
- The system is suitable also for other anaerobic studies such as Denitrification Tests (Requires 6 Denitrification glass bottles A00000410).
- Easy-to-use and powerful sensor with 4 digit LED Display.
- Unique wireless control from your PC and VELP Ermes cloud platform.



Wireless DataBox™
and RESPIROSoft™

INSTRUMENT - CODE

| | | |
|--|------------------|------------|
| RESPIROMETRIC Sensor System 6 Maxi | 230 V / 50-60 Hz | SA10200156 |
| RESPIROMETRIC Sensor System 6 Maxi | 115 V / 60 Hz | SA10210156 |
| RESPIROMETRIC Sensor System 6 Maxi no Databox™ | 230 V / 50-60 Hz | SB10200156 |
| RESPIROMETRIC Sensor System 6 Maxi no Databox™ | 115 V / 60 Hz | SB10210156 |

Included in the system: Maxi stirring station, 1000 ml bottles, KOH containers, stirring bars, DataBox™ and RESPIROSoft™ software.

ermes enabled



RESPIROMETRIC Sensor System 6 for Soil Analysis

The RESPIROMETRIC Sensor System for soil analysis is designed to support aerobic degradation tests on soil and compost. The analysis provides microbial soil respiration results from the mineralization of organic substances. A major application of these test is related to the detection and monitoring processes during decontaminations and clean-ups.

- The analysis is carried out directly on soil/compost in the 1000ml flask.
- Results reading in Pressure values (millibar).
- Easy-to-use and powerful sensor with 4 digit LED Display.
- Unique wireless control from your PC and VELP Ermes cloud platform.



INSTRUMENT - CODE

| | | |
|---|---|------------|
| RESPIROMETRIC Sensor System for Soil Analysis | - | SA102B0176 |
| RESPIROMETRIC Sensor System for Soil Analysis no Databox™ | - | SB102B0176 |

Included in the system: 6 Kits soil analysis, DataBox™ and RESPIROSoft™ software.



RESPIROMETRIC Sensor System 6 for Plastic Biodegradability

The RESPIROMETRIC Sensor System for plastic biodegradability is designed to support aerobic degradation test for the determination of the degree of aerobic biodegradability of plastic materials, including those containing formulation additives, in an aqueous medium.

- Multistirrer Digital 6-positions to precisely set stirring speed and enable the lock feature to avoid unintentional modification of the stirring speed.
- Results reading in Pressure values (millibar).
- Easy-to-use and powerful sensor with 4 digit LED Display.
- Unique wireless control from your PC and VELP Ermes cloud platform.



INSTRUMENT - CODE

| | | |
|--|---------------------|------------|
| RESPIROMETRIC Sensor System for Plastic Biodegradability | 100-240 V / 50-60Hz | SA102A0166 |
| RESPIROMETRIC Sensor System for Plastic Biodegradability no Databox™ | 100-240 V / 50-60Hz | SB102A0166 |

Included in the system: Multistirrer Digital 6-positions, 6 Kits plastic biodegradability analysis, stirring bars, DataBox™ and RESPIROSoft™ software.



■ BOD Sensor

Innovative, mercury-free and extremely reliable solution for BOD analysis. The 6 and 10 position systems include the BOD sensors, stirring station, stirring bars and KOH traps.

PRECISE AND EASY TO USE

It can measure on scales of 90, 250, 600 and 999 ppm BOD.

SAFE

The IP 54 Protection Rating guarantee durability even in the most aggressive environments. Class 3 IEC 1010 conformity.

COMPACT DESIGN

BOD System 6 is designed to optimize space consumption inside cooled incubators.



INSTRUMENT - CODE

| | | |
|----------------|---|-----------|
| BOD Sensor | - | F102B0133 |
| BOD Sensor Set | - | F102B0134 |



■ BOD Sensor System 6

Innovative, mercury-free and extremely reliable solution for BOD analysis. The 6 and 10 position systems include the BOD sensors, stirring station, stirring bars and KOH traps.

PRECISE AND EASY TO USE

It can measure on scales of 90, 250, 600, 999 and 4000 ppm BOD.

SAFE

The IP 54 Protection Rating guarantee durability even in the most aggressive environments. Class 3 IEC 1010 conformity.

COMPACT DESIGN

BOD Sensor System 6 is designed to optimize space consumption inside cooled incubators.



BOD Scales
90, 250, 600, 999, 4000 ppm

INSTRUMENT - CODE

| | | |
|---------------------|------------------|------------|
| BOD Sensor System 6 | 230 V / 50-60 Hz | SA10200136 |
| BOD Sensor System 6 | 115 V / 60 Hz | SA10210136 |

Included in the system: Stirring station 6 positions, BOD Sensor, 500 ml amber bottles, KOH containers, stirring bars.



BOD Sensor System 10

The VELP RESPIROMETRIC Sensor is designed to provide a smart and cloud-enabled solution to laboratories of water and wastewater quality control, biological waste treatment plants, renewable energy production and environmental labs studying the effects of plastic and other pollutants in water and soils.

PRECISE AND CONVENIENT

It can measure on scales of 90, 250, 600, 999 and 4000 ppm BOD.

COMPACT DESIGN

BOD Sensor System 10 is designed to optimize space consumption inside cooled incubators.



BOD Scales
90, 250, 600, 999, 4000 ppm

INSTRUMENT - CODE

| | | |
|----------------------|------------------|------------|
| BOD Sensor System 10 | 230 V / 50-60 Hz | SA10200137 |
| BOD Sensor System 10 | 115 V / 60 Hz | SA10210137 |

Included in the system: Stirring station 10 positions, BOD Sensor, 500 ml amber bottles, KOH containers, stirring bars.



ACCESSORIES - BOD AND RESPIROMETERS

| | RESPIROMETRIC Sensor System 6 BOD | RESPIROMETRIC Sensor System 6 Maxi BMP | RESPIROMETRIC Sensor System 6 for Soil Analysis | RESPIROMETRIC Sensor System 6 for Plastic Biodegradability | BOD Sensor System |
|---|-----------------------------------|--|---|--|-------------------|
| A00000135 Sensor Check | • | • | • | • | • |
| A00000136 Control Test Tablets, 10 pcs/box | • | • | • | • | • |
| A00000221 Multi-socket extension cable EU | • | • | | • | • |
| A00000222 Multi-socket extension cable UK | • | • | | • | • |
| A00000223 Multi-socket extension cable AU | • | • | | • | • |
| 10003083 Plug USA for power supply | • | • | | • | • |
| 10003084 Plug UK for power supply | • | • | | • | • |
| 10003085 Plug AU for power supply | • | • | | • | • |
| A00000410 Denitrification glass bottle | | • | | | |
| E00010012 VELP Ermes 1 year Connection | • | • | • | • | |
| E00010036 VELP Ermes 3 year Connection | • | • | • | • | |

Cooled Incubators

■ FOC 200IL Connect (not available in USA, Mexico and Canada)

FOC 200IL Connect is the cloud-enabled illuminated cooled incubator ideal for plant growth, seed germination, cultivation, breeding of plants, insects and small animals, bacteria culturing, fermentation studies.

CONVENIENT AND EFFICIENT

FOC 200IL Connect is equipped with 2 lighted shelves and an internal transparent door.

EFFICIENT THERMOREGULATION

FOC 200IL Connect guarantees a perfect thermoregulation and stability thanks to the Auto-Tuning system. The use of forced ventilation uniforms the internal temperature.

EXCLUSIVE CONNECTIVITY

Premium Wi-Fi connectivity to VELP Ermes Cloud Platform. TEMPSoft™ software with 485 cable connection to constant monitor the instrument. MODBUS communications for PLC connection.



ermes enabled



Total volume
169 L



Temperature range
From 3 to 50 °C



Auto-Tuning
technology
± 0.5°C



Timer
Lighting Cycles

INSTRUMENT - CODE

FOC 200IL Connect

230 V / 50-60 Hz

F10400544*

* The cooled incubator is supplied with 2 shelves 40002722 and 1 shelf 10008039



■ FOC 200I Connect (not available in USA, Mexico and Canada)

FOC 200I Connect is the cloud-enabled cooled incubator with internal transparent door ensuring excellent temperature stability and homogeneity.

CONVENIENT AND EFFICIENT

FOC 200I Connect is equipped with an internal transparent door allowing sample monitoring with no temperature alteration.

EFFICIENT THERMOREGULATION

FOC 200I Connect guarantees a perfect thermoregulation and stability thanks to the Auto-Tuning system. The use of forced ventilation uniforms the internal temperature.

EXCLUSIVE CONNECTIVITY

Premium Wi-Fi connectivity to VELP Ermes Cloud Platform. TEMPSoft™ software with 485 cable connection to constant monitor the instrument. MODBUS communications for PLC connection.



ermes enabled



Total volume
169 L



Temperature range
From 3 to 50 °C



Auto-Tuning
technology
± 0.5°C

INSTRUMENT - CODE

FOC 200I Connect

230 V / 50-60 Hz

F10400543*

* The cooled incubator is supplied with 3 shelves 10008039



■ FOC 200E Connect (not available in USA, Mexico and Canada)

FOC 200E Connect is the cloud-enabled cooled incubator ensuring excellent temperature stability and homogeneity.

CONVENIENT AND EFFICIENT

FOC 200E Connect is easy to use thanks to the intuitive digital interface.

EFFICIENT THERMOREGULATION

FOC 200E Connect guarantees a perfect thermoregulation and stability thanks to the Auto-Tuning system. The use of forced ventilation uniforms the internal temperature.

EXCLUSIVE CONNECTIVITY

Premium Wi-Fi connectivity to VELP Ermes Cloud Platform. TEMPSoft™ software with 485 cable connection to constantly monitor the instrument. MODBUS communications for PLC connection.



ermes enabled



Total volume
169 L



Temperature range
From 3 to 50 °C



Auto-Tuning
technology
± 0.5°C

INSTRUMENT - CODE

FOC 200E Connect

230 V / 50-60 Hz

F10300542*

* The cooled incubator is supplied with 3 shelves 10008039



■ FOC 120I Connect (not available in USA, Mexico and Canada)

FOC 120I Connect is the cloud-enabled cooled incubator with internal transparent door ensuring excellent temperature stability and homogeneity.

CONVENIENT AND EFFICIENT

FOC 120I Connect is equipped with an internal transparent door allowing sample monitoring with no temperature alteration.

EFFICIENT THERMOREGULATION

FOC 120I Connect guarantees a perfect thermoregulation and stability thanks to the Auto-Tuning system. The use of forced ventilation uniforms the internal temperature.

EXCLUSIVE CONNECTIVITY

Premium Wi-Fi connectivity to VELP Ermes Cloud Platform. TEMPSoft™ software with 485 cable connection to constantly monitor the instrument. MODBUS communications for PLC connection.



ermes enabled



Total volume
109 L



Temperature range
From 3 to 50 °C



Auto-Tuning
technology
± 0.5°C

INSTRUMENT - CODE

FOC 120I Connect

230 V / 50-60 Hz

F10400541

* The cooled incubator is supplied with 2 shelves



■ FOC 120E Connect (not available in USA, Mexico and Canada)

FOC 120E Connect is the cloud-enabled cooled incubator ensuring excellent temperature stability and homogeneity.

CONVENIENT AND EFFICIENT

FOC 120E Connect is asy to use thanks to the intuitive digital interface.

EFFICIENT THERMOREGULATION

FOC 120E Connect guarantees a perfect thermoregulation and stability thanks to the Auto-Tuning system. The use of forced ventilation uniforms the internal temperature.

EXCLUSIVE CONNECTIVITY

Premium Wi-Fi connectivity to VELP Ermes Cloud Platform. TEMPSoft™ software with 485 cable connection to constant monitor the instrument. MODBUS communications for PLC connection.



ermes enabled



Total volume
109 L



Temperature range
From 3 to 50 °C



Auto-Tuning
technology
± 0.5°C

INSTRUMENT - CODE

FOC 120E Connect

230 V / 50-60 Hz

F10300540

* The cooled incubator is supplied with 2 shelves



■ FTC 120 (not available in USA, Mexico and Canada)

Compact and efficient incubator specially designed for BOD respirometric measurement ensuring a constant 20°C temperature. FTC 120 is equipped with 2 internal sockets.

EFFICIENT THERMOREGULATION

The use of forced ventilation uniforms the internal temperature.



Total volume
109 L



Temperature
20 °C



Auto-Tuning
technology
± 0.5°C



INSTRUMENT - CODE

FTC 120

230 V / 50-60 Hz

F10300305

* The cooled incubator is supplied with 2 shelves



ACCESSORIES - COOLED INCUBATORS

| | FOC 200IL Connect | FOC 200I Connect | FOC 200E Connect | FOC 120I Connect | FOC 120E Connect | FTC 120 |
|---|----------------------|---------------------|---------------------|---------------------|---------------------|---------|
| A00000404 TEMPSoft with RS485 cable | • | • | • | • | • | |
| A00000405 Interconnection cable FOC Connect | • | • | • | • | • | |
| A00000406 IQ/OQ FOC Connect and FTC 120 Manual | • | • | • | • | • | • |
| A00000221 Multi-socket extension cable EU | • | • | • | • | • | • |
| A00000222 Multi-socket extension cable UK | • | • | • | • | • | • |
| A00000223 Multi-socket extension cable AU | • | • | • | • | • | • |
| 10008039 Shelf for FOC 200, FOC 120, FTC 120 | • | • | • | • | • | • |
| 10008040 Bottom Shelf FOC 120, FTC 120 | | | | • | • | • |
| 40002722 Lighted shelf for FOC 200 IL | • | | | | | |
| E00010012 VELP Ermes 1 year Connection | • | • | • | • | • | |
| E00010036 VELP Ermes 3 year Connection | • | • | • | • | • | |



Flocculators

■ JLT 6

Easy to use and reliable digital flocculators with 6 positions, timer and illuminated back panel. JLT 6 is the ideal solution for Jar Testing 6 samples with the same speed and time settings.

RELIABLE

JLT guarantees a precise speed thanks to the electronically controlled motor. It also ensures a strong resistance to chemical and mechanical corrosion.

PRACTICAL

JLT has an illuminated back panel that simplifies the observation of the samples. All the stirring rods are adjustable in height without tools.



Stirring speed
From 10 to 300 rpm



Timer
From 0 to 999 min
From 0 to 99 hours



INSTRUMENT - CODE

JLT 6

100-240 V / 50-60 Hz

F105A0109



■ JLT 4

Easy to use and reliable digital flocculators with 4 positions, timer and illuminated back panel. JLT 4 is the ideal solution for Jar Testing 4 samples with the same speed and time settings.

RELIABLE

JLT guarantees a precise speed thanks to the electronically controlled motor. It also ensures a strong resistance to chemical and mechanical corrosion.

PRACTICAL

JLT has an illuminated back panel that simplifies the observation of the samples. All the stirring rods are adjustable in height without tools.



Stirring speed
From 10 to 300 rpm



Timer
From 0 to 999 min
From 0 to 99 hours



INSTRUMENT - CODE

JLT 4

100-240 V / 50-60 Hz

F105A0108



■ FC6S

Precise and easy to use flocculators with 6 independent positions that allow standard conditions to be adopted for Jar Test. FC6S has an illuminated back panel and adjustable stirring rods.

RESISTANT

Strong resistance to chemical and mechanical corrosion. Precise speed setting with the electronically controlled motor.

VERSATILE

With the FC6S, the same product can be submitted to different conditions in parallel, while the speed can be changed during the run.



Stirring speed steps
10 - 15 - 30 - 45 - 60 - 90
120 - 150 - 200 - 300 rpm



INSTRUMENT - CODE

FC6S

100-240 V / 50-60 Hz

F105A0112



■ FC4S

Precise and easy to use flocculators with 4 independent positions that allow standard conditions to be adopted for Jar Test. FC4S has an illuminated back panel and adjustable stirring rods.

RESISTANT

Strong resistance to chemical and mechanical corrosion. Precise speed setting with the electronically controlled motor.

VERSATILE

With the FC4S, the same product can be submitted to different conditions in parallel, while the speed can be changed during the run.



Stirring speed steps
10 - 15 - 30 - 45 - 60 - 90
120 - 150 - 200 - 300 rpm



INSTRUMENT - CODE

FC4S

100-240 V / 50-60 Hz

F105A0111



■ FP4

Portable flocculator with anti-skid base for stability, that allows standard conditions to be adopted for Jar Test. FP4 allows 4 samples to be processed in parallel. FP4 can be connected to a battery or a car cigarette lighter (12V).

RELIABLE

FP4 guarantees a precise speed thanks to an electronically controlled DC gear motor. It also ensures a strong resistance to chemical and mechanical corrosion.

PRACTICAL

FP4 is designed to be used on-site for immediate results. All the stirring rods are adjustable in height without tools.



Stirring speed steps
20 - 40 - 50 - 100 - 200 rpm



Timer
From 0 to 30 min

INSTRUMENT - CODE

FP 4

100-240 V / 50-60 Hz

F105A0117



ACCESSORIES - FLOCCULATORS

| | JLT 6 | JLT 4 | FC6S | FC4S | FP4 |
|--|-------|-------|------|------|-----|
| A00001000 Plastic beaker 1000 ml | • | • | • | • | • |
| A00001001 Glass beaker 1000 ml | • | • | • | • | • |
| A00001002 Transparent plastic Imhoff cone | • | • | • | • | |
| A00001003 Glass graduated Imhoff cone | • | • | • | • | |
| A00001004 Stand for 2 Imhoff cones | • | • | • | • | |
| A00001005 Carrying case | | | | | • |
| A00000319 IQ/OQ Manual JLT | • | • | | | |
| A00000379 IQ/OQ FC4S - FC6S Manual | | | • | • | |
| 10003722 Power supply cord, US plug | • | • | • | • | |
| 10003723 Power supply cord, UK plug | • | • | • | • | |
| 10003724 Power supply cord, AU plug | • | • | • | • | |
| 10003083 Plug USA for power supply | | | | | • |
| 10003084 Plug UK for power supply | | | | | • |
| 10003085 Plug AU for power supply | | | | | • |

Other Solutions

■ ROTAX 6.8 Overhead Mixer

Overhead mixer with LCD speed display, that offers excellent performance in applications where rotating, shaking and mixing are required. Rotax 6.8 is resistant to chemical and mechanical corrosion.

ACCURATE AND SAFE

Rotax 6.8 is designed to automatically stop if blocked. Moreover it is compliant with UNI 10802 and DIN 38414 standards.



Stirring speed
Up to 30 rpm



INSTRUMENT - CODE

ROTAX 6.8

230 V / 50-60 Hz

F10600118

ROTAX 6.8

115 V / 50-60 Hz

F10610118

For accessories visit online www.velp.com



■ OCB Open Circulating Baths

Efficient and precise open circulating bath that ensures excellent temperature homogeneity ($\pm 1^\circ\text{C}$), accurate control and thermostating of testing materials.

EASY TO USE

LCD display ensures a precise temperature setting. A wide choice of test tube racks is available.



Temperature range
From Room Temperature to 105 °C



Timer
From 0 to 24:59 hours



INSTRUMENT - CODE

OCB

230 V / 50-60 Hz

F40300240

OCB

115 V / 50-60 Hz

F40310240

For accessories visit online www.velp.com



■ Radiation Detector (not available in USA, Mexico and Canada)

Standard model of a Geiger-Müller tube with mica output window, able to detect α , β , γ and x-rays. Results are shown by a graduated scale, flashing LED and acoustic signal

EASY TO USE AND MULTIPURPOSE

It ensures an operating time of up to 2000 hours with one set of batteries and gives results in milliroentgens/hour and counts/minute (cpm) in 3 periods.

RELIABLE

The anti-saturation circuit reaches a value equivalent to 100 times the maximum reading in the maximum range.



INSTRUMENT - CODE

Radiation detector

For accessories visit online www.velp.com

R10800340



■ SP 311 Peristaltic Pumps

An easy to use and reliable peristaltic pump, ideal for transferring gas, liquids, mixtures and emulsions in various industries, able to pump from 10 meters depth.

EASY TO USE AND RELIABLE

Tubing resists most chemicals and ensures minimum maintenance. The electronic speed control makes simple the control of the instrument.



INSTRUMENT - CODE

| | | |
|-----------|------------------|-----------|
| SP 311/2 | 230 V / 50-60 Hz | F40200002 |
| SP 311/6 | 230 V / 50-60 Hz | F40200006 |
| SP 311/12 | 230 V / 50-60 Hz | F40200012 |
| SP 311/60 | 230 V / 50-60 Hz | F40200060 |

For accessories visit online www.velp.com



■ TB1 Turbidimeter (not available in USA, Mexico and Canada)

Waterproof, precise and easy to calibrate portable turbidimeter that measures the turbidity of aqueous samples simply and accurately with nephelometric method.

RELIABLE AND PRECISE

Tb1 is supplied with 4 calibration standards (800, 100, 20 and 0.02 NTU) and can test up to 1200 samples with one set of batteries.

RESISTANT

IP 67 protection against ingress of particles and water.



INSTRUMENT - CODE

TB 1

For accessories visit online www.velp.com

R109B12150



■ JPV Recirculating Water Vacuum Pump

JPV pump ensures constant water recirculation and reduced water consumption. JPV resists to chemical corrosion and is ideal for applications where non-pulsing suction is needed.

QUALITY AND USABILITY

JPV provides vacuum where no reliable water supply is available and guarantees the absence of noise and vibrations. ABS structure guarantees high resistance to chemical corrosion.

EASY TO USE

The instrument has inset handles to simplify lifting and the stopcock for emptying the tank.



INSTRUMENT - CODE

JPV

230 V / 50-60 Hz

F40100541

JPV

115 V / 50-60 Hz

F40110541

For accessories visit online www.velp.com





ITALY – HQ

Via Stazione 16 20865 Usmate (MB) Italy
Tel. +39 039 628811 - velpitalia@velp.com

USA

40, Burt Drive, Unit #1, Deer Park - NY 11729 - U.S.
Tel. +1 631 573 6002 - velpusa@velp.com

CHINA

Building 28, Lane 1333, Xinlong road, Minhang district, Shanghai
Tel. +8621 34500630 - velpchina@velp.com

INDIA

velpindia@velp.com

LATAM

velplatam@velp.com



Rev.1.5.2.2022

VELP Official Partner

We reserve the right to make technical alterations
We do not assume liability for errors in printing, typing or transmission



DESIGNED AND MANUFACTURED IN ITALY